



EBERLINE SERVICES

EBERLINE ANALYTICAL CORPORATION
2030 Wright Avenue
Richmond, California 94804-3849
Phone (510) 235-2633 Fax (510) 235-0438
Toll Free (800) 841-5487
www.eberlineservices.com

June 13, 2008

Mr. Steve Trent
Fluor Hanford Inc.
1200 Jadwin Avenue
Richland, WA 99352



Reference: **P.O. #33677**
Eberline Services R8-04-167-7087, SDG H3702

Dear Mr. Trent:

Enclosed is the data report for one solid (soil) sample designated under SAF No. F08-070 received at Eberline Services on April 29, 2008. The sample was analyzed according to the accompanying chain-of-custody document.

Please call if you have any questions concerning this report.

Sincerely,

Melissa C. Mannion
Senior Program Manager

MCM/njv

Enclosure: Data Package

00000001

1.0 GENERAL

Fluor Hanford Inc. (FH) Sample Delivery Group H3702 was composed of one solid (soil) sample designated under SAF No. F08-070 with a Project Designation of: 200-TW-1 OU Characterization for Well 299-E33-342 – Soil.

The sample was received as stated on the Chain-of-Custody document. Any discrepancies are noted on the Eberline Services Sample Receipt Checklist.

2.0 ANALYSIS NOTES

2.1 Tritium Analysis

No problems were encountered during the course of the analyses.

2.2 Carbon-14 Analysis

No problems were encountered during the course of the analyses.

2.3 Nickel-63 Analysis

No problems were encountered during the course of the analyses.

2.4 Technetium-99 Analysis

No problems were encountered during the course of the analyses.

2.5 Iodine-129 Analysis

No problems were encountered during the course of the analyses.

2.6 Radium-226 Analysis

No problems were encountered during the course of the analyses.

2.7 Radium-228 Analysis

No problems were encountered during the course of the analyses.

2.8 Isotopic Thorium Analysis


No problems were encountered during the course of the analyses.

3.0 Case Narrative Certification Statement

"I certify that this data package is in compliance with the SOW, both technically and for completeness, for other than the conditions detailed above. Release of the data obtained in this hard copy data package has been authorized by the Laboratory Manager or a designee, as verified by the following signature."



Melissa C. Mannion
Senior Program Manager

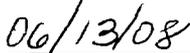

Date

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

Client Hanford
Contract No. 33677
Case no SDG_H3702

S U M M A R Y D A T A S E C T I O N

T A B L E O F C O N T E N T S				
About this section	.	.	.	1
Sample Summaries	.	.	.	3
Prep Batch Summary	.	.	.	5
Work Summary	.	.	.	6
Method Blanks	.	.	.	8
Lab Control Samples	.	.	.	9
Duplicates	.	.	.	10
Data Sheets	.	.	.	11
Method Summaries	.	.	.	12
Report Guides	.	.	.	20
End of Section	.	.	.	34

Prepared by



Reviewed by



Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-TOC
Version 3.06
Report date 06/13/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG H3702

ABOUT THE DATA SUMMARY SECTION

The Data Summary Section of a Data Package has all data, in several useful orders, necessary for first level, routine review of the data package for a Sample Delivery Group (SDG). This section follows the Data Package Narrative, which has an overview of the data package and a discussion of special problems. It is followed by the Raw Data Section, which has full details.

The Data Summary Section has several groups of reports:

SAMPLE SUMMARIES

The Sample and QC Summary Reports show all samples, including QC samples, reported in one SDG. These reports cross-reference client and lab sample identifiers.

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches (lab groupings reflecting how work was organized) relevant to the reported SDG with information necessary to check the completeness and consistency of the SDG.

WORK SUMMARY

The Work Summary Report shows all samples and work done on them relevant to the reported SDG.

METHOD BLANKS

The Method Blank Reports, one for each Method Blank relevant to the SDG, show all results and primary supporting information for the blanks.

LAB CONTROL SAMPLES

The Lab Control Sample Reports, one for each Lab Control Sample relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 1

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000005

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. 33677
Case no SDG_H3702

ABOUT THE DATA SUMMARY SECTION

DUPLICATES

The Duplicate Reports, one for each Duplicate and Original sample pair relevant to the SDG, show all results, differences and primary supporting information for these QC samples.

MATRIX SPIKES

The Matrix Spike Reports, one for each Spiked and Original sample pair relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

DATA SHEETS

The Data Sheet Reports, one for each client sample in the SDG, show all results and primary supporting information for these samples.

METHOD SUMMARIES

The Method Summary Reports, one for each test used in the SDG, show all results, QC and method performance data for one analyte on one or two pages. (A test is a short code for the method used to do certain work to the client's specification.)

REPORT GUIDES

The Report Guides, one for each of the above groups of reports, have documentation on how to read the associated reports.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 2

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000006

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

LAB SAMPLE SUMMARY

Client Hanford
 Contract No. 33677
 Case no SDG H3702

LAB SAMPLE ID	CLIENT SAMPLE ID	LOCATION	MATRIX	LEVEL	SAF NO	CHAIN OF CUSTODY	COLLECTED
R804167-01	B1TNRO	C5857, I-098	SOLID		F08-070	F08-070-028	04/01/08 08:20
R804167-02	Lab Control Sample		SOLID		F08-070		
R804167-03	Method Blank		SOLID		F08-070		
R804167-04	Duplicate (R804167-01)	C5857, I-098	SOLID		F08-070		04/01/08 08:20

LAB SUMMARY

Page 1

SUMMARY DATA SECTION

Page 3

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LS
 Version 3.06
 Report date 06/13/08

00000007

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

QC SUMMARY

Client Hanford
 Contract No. 33677
 Case no SDG H3702

QC BATCH	CHAIN OF CUSTODY	CLIENT SAMPLE ID	MATRIX	% SOLIDS	SAMPLE AMOUNT	BASIS AMOUNT	DAYS SINCE RECEIVED	LAB COLL	DEPARTMENT SAMPLE ID
7087	F08-070-028	B1TNR0	SOLID	89.3	157 g		04/29/08 28	R804167-01	7087-001
		Method Blank	SOLID					R804167-03	7087-003
		Lab Control Sample	SOLID					R804167-02	7087-002
		Duplicate (R804167-01)	SOLID	89.3	157 g		04/29/08 28	R804167-04	7087-004

QC SUMMARY

Page 1

SUMMARY DATA SECTION

Page 4

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-QS
 Version 3.06
 Report date 06/13/08

00000008

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

PREP BATCH SUMMARY

Client Hanford
 Contract No. 33677
 Case no SDG H3702

TEST	MATRIX	METHOD	PREPARATION ERROR		PLANCHETS ANALYZED				QUALI-	
			BATCH	2σ %	CLIENT	MORE	RE	BLANK		LCS
Alpha Spectroscopy										
TH	SOLID	Thorium, Isotopic in Solids	6150-050	8.0	1			1	1	1/1
Beta Counting										
AC	SOLID	Radium 228 in Solids	6150-050	10.4	1			1	1	1/1
TC	SOLID	Technetium 99 in Solids	6150-050	13.2	1			1	1	1/1
Gamma Spectroscopy										
I	SOLID	Iodine 129 in Solids	6150-050	19.4	1			1	1	1/1
Liquid Scintillation Counting										
C	SOLID	Carbon 14 in Solids	6150-050	10.0	1			1	1	1/1
H	SOLID	Tritium in Solids	6150-050	10.0	1			1	1	1/1
NI_L	SOLID	Nickel 63 in Solids	6150-050	11.2	1			1	1	1/1
Radon Counting										
RA	SOLID	Radium 226 in Solids	6150-050	16.4	1			1	1	1/1

Duplicates and Matrix Spikes are those with original (Client) sample in this Sample Delivery Group.
 Blank and LCS planchets are those in the same preparation batch as some Client, Duplicate or Spike sample.

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-PBS
 Version 3.06
 Report date 06/13/08

00000009

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

LAB WORK SUMMARY

Client Hanford
 Contract No. 33677
 Case no SDG H3702

LAB SAMPLE	CLIENT SAMPLE ID				SUF-					
COLLECTED	LOCATION	MATRIX			FIX	ANALYZED	REVIEWED	BY	METHOD	
RECEIVED	CUSTODY	SAF No	PLANCHET	TEST						
R804167-01	B1TNR0		7087-001	AC		06/09/08	06/13/08	BW	Radium 228 in Solids	
04/01/08	C5857, I-098	SOLID	7087-001	C		06/02/08	06/05/08	BW	Carbon 14 in Solids	
04/29/08	F08-070-028	F08-070	7087-001	H		05/30/08	06/03/08	BW	Tritium in Solids	
			7087-001	I		06/05/08	06/09/08	BW	Iodine 129 in Solids	
			7087-001	NI_L		05/28/08	06/03/08	BW	Nickel 63 in Solids	
			7087-001	RA		05/27/08	05/29/08	BW	Radium 226 in Solids	
			7087-001	TC		05/29/08	06/13/08	BW	Technetium 99 in Solids	
			7087-001	TH		06/03/08	06/06/08	BW	Thorium, Isotopic in Solids	
R804167-02	Lab Control Sample		7087-002	AC		06/09/08	06/13/08	BW	Radium 228 in Solids	
		SOLID	7087-002	C		06/02/08	06/05/08	BW	Carbon 14 in Solids	
		F08-070	7087-002	H		05/30/08	06/03/08	BW	Tritium in Solids	
			7087-002	I		06/05/08	06/09/08	BW	Iodine 129 in Solids	
			7087-002	NI_L		05/28/08	06/03/08	BW	Nickel 63 in Solids	
			7087-002	RA		05/27/08	05/29/08	BW	Radium 226 in Solids	
			7087-002	TC		05/27/08	06/13/08	BW	Technetium 99 in Solids	
			7087-002	TH		06/05/08	06/06/08	BW	Thorium, Isotopic in Solids	
R804167-03	Method Blank		7087-003	AC		06/09/08	06/13/08	BW	Radium 228 in Solids	
		SOLID	7087-003	C		06/02/08	06/05/08	BW	Carbon 14 in Solids	
		F08-070	7087-003	H		05/30/08	06/03/08	BW	Tritium in Solids	
			7087-003	I		06/06/08	06/09/08	BW	Iodine 129 in Solids	
			7087-003	NI_L		05/28/08	06/03/08	BW	Nickel 63 in Solids	
			7087-003	RA		05/27/08	05/29/08	BW	Radium 226 in Solids	
			7087-003	TC		05/27/08	06/13/08	BW	Technetium 99 in Solids	
			7087-003	TH		06/05/08	06/06/08	BW	Thorium, Isotopic in Solids	
R804167-04	Duplicate (R804167-01)		7087-004	AC		06/09/08	06/13/08	BW	Radium 228 in Solids	
04/01/08	C5857, I-098	SOLID	7087-004	C		06/02/08	06/05/08	BW	Carbon 14 in Solids	
04/29/08		F08-070	7087-004	H		05/30/08	06/03/08	BW	Tritium in Solids	
			7087-004	I		06/06/08	06/09/08	BW	Iodine 129 in Solids	
			7087-004	NI_L		05/28/08	06/03/08	BW	Nickel 63 in Solids	
			7087-004	RA		05/27/08	05/29/08	BW	Radium 226 in Solids	
			7087-004	TC		05/29/08	06/13/08	BW	Technetium 99 in Solids	
			7087-004	TH		06/03/08	06/06/08	BW	Thorium, Isotopic in Solids	

WORK SUMMARY

Page 1

SUMMARY DATA SECTION

Page 6

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LWS
 Version 3.06
 Report date 06/13/08

00000010

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087

Contact Melissa C. Mannion

WORK SUMMARY, cont.

Client Hanford

Contract No. 33677

Case no SDG H3702

COUNTS OF TESTS BY SAMPLE TYPE

TEST	SAF No	METHOD	REFERENCE	CLIENT	MORE	RE	BLANK	LCS	DUP SPIKE	TOTAL
AC	F08-070	Radium 228 in Solids	RAISO_SEP_GPC	1			1	1	1	4
C	F08-070	Carbon 14 in Solids	C14_COX_LSC	1			1	1	1	4
H	F08-070	Tritium in Solids	TRITIUM_COX_LSC	1			1	1	1	4
I	F08-070	Iodine 129 in Solids	I129_SEP_LEPS_GS	1			1	1	1	4
NI_L	F08-070	Nickel 63 in Solids	NI63_LSC	1			1	1	1	4
RA	F08-070	Radium 226 in Solids	903.1_RA226_LUC	1			1	1	1	4
TC	F08-070	Technetium 99 in Solids	TC99_TR_SEP_GPC	1			1	1	1	4
TH	F08-070	Thorium, Isotopic in Solids	THISO_IE_PLATE_AEA	1			1	1	1	4
TOTALS				8			8	8	8	32

WORK SUMMARY

Page 2

SUMMARY DATA SECTION

Page 7

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LWS
 Version 3.06
 Report date 06/13/08

00000011

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP H3702

7087-003

Method Blank

METHOD BLANK

SDG <u>7087</u>	Client/Case no <u>Hanford</u>	SDG <u>H3702</u>
Contact <u>Melissa C. Mannion</u>	Contract <u>No. 33677</u>	
Lab sample id <u>R804167-03</u>	Client sample id <u>Method Blank</u>	
Dept sample id <u>7087-003</u>	Material/Matrix <u>SOLID</u>	
	SAF No <u>F08-070</u>	

ANALYTE	CAS NO	RESULT pCi/g	2σ ERR (COUNT)	MDA pCi/g	RDL pCi/g	QUALI- FIERS	TEST
Tritium	10028-17-8	0.514	2.8	4.78	400	U	H
Carbon 14	14762-75-5	2.12	2.3	3.75	50.0	U	C
Nickel 63	13981-37-8	-1.45	1.8	3.16	30.0	U	NI_L
Technetium 99	14133-76-7	-0.005	0.14	0.486	12.0	U	TC
Thorium 228	14274-82-9	0.156	0.25	0.418	1.00	U	TH
Thorium 230	14269-63-7	0.031	0.19	0.237	1.00	U	TH
Thorium 232	TH-232	-0.031	0.062	0.237	1.00	U	TH
Iodine 129	15046-84-1	-0.273	0.58	1.32	2.00	U	I
Radium 226	13982-63-3	0	0.037	0.071	0.100	U	RA
Radium 228	15262-20-1	-0.192	0.35	0.978	1.00	U	AC

200TW1 OU Charactrzn 299E33342-Soil

QC-BLANK #65614

Lab id <u>EBRLNE</u>
Protocol <u>Fluor</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>06/13/08</u>

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

7087-002

Lab Control Sample

LAB CONTROL SAMPLE

SDG <u>7087</u> Contact <u>Melissa C. Mannion</u> Lab sample id <u>R804167-02</u> Dept sample id <u>7087-002</u>	Client/Case no <u>Hanford</u> SDG <u>H3702</u> Contract No. <u>33677</u> Client sample id <u>Lab Control Sample</u> Material/Matrix _____ <u>SOLID</u> SAF No <u>F08-070</u>
---	--

ANALYTE	RESULT	2σ ERR	MDA	RDL	QUALI-	ADDED	2σ ERR	REC	3σ LMTS	PROTOCOL
	pCi/g	(COUNT)	pCi/g	pCi/g	FIERS TEST	pCi/g	pCi/g	%	(TOTAL)	LIMITS
Tritium	581	14	5.25	400	H	615	25	94	84-116	80-120
Carbon 14	1480	14	3.46	50.0	C	1600	64	92	85-115	80-120
Nickel 63	246	6.1	3.15	30.0	NI_L	266	11	92	83-117	80-120
Technetium 99	110	2.4	0.542	12.0	TC	120	4.8	92	81-119	80-120
Thorium 230	20.7	2.2	0.214	1.00	TH	21.9	0.88	95	80-120	80-120
Iodine 129	134	1.5	1.88	2.00	I	127	5.1	106	69-131	80-120
Radium 226	6.41	0.32	0.082	0.100	RA	5.58	0.22	115	70-130	80-120
Radium 228	8.88	0.66	0.967	1.00	AC	9.52	0.38	93	81-119	80-120

200TW1 OU Charactrzn 299E33342-Soil

QC-LCS #65613

Lab id <u>EBRLNE</u>
Protocol <u>Fluor</u>
Version <u>Ver 1.0</u>
Form <u>DVD-LCS</u>
Version <u>3.06</u>
Report date <u>06/13/08</u>

00000013

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

7087-004

B1TNR0

DUPLICATE

SDG <u>7087</u>	Client/Case no <u>Hanford</u>	SDG <u>H3702</u>
Contact <u>Melissa C. Mannion</u>	Contract No. <u>33677</u>	
DUPLICATE	ORIGINAL	
Lab sample id <u>R804167-04</u>	Lab sample id <u>R804167-01</u>	Client sample id <u>B1TNR0</u>
Dept sample id <u>7087-004</u>	Dept sample id <u>7087-001</u>	Location/Matrix <u>C5857, I-098</u> <u>SOLID</u>
	Received <u>04/29/08</u>	Collected/Weight <u>04/01/08 08:20</u> <u>157 g</u>
% solids <u>89.3</u>	% solids <u>89.3</u>	Custody/SAF No <u>F08-070-028</u> <u>F08-070</u>

ANALYTE	DUPLICATE	2σ ERR	MDA	RDL	QUALI-	ORIGINAL	2σ ERR	MDA	QUALI-	RPD	3σ	DER
	pCi/g	(COUNT)	pCi/g	pCi/g	FIERS TEST		pCi/g	(COUNT)	pCi/g	FIERS	%	TOT
Tritium	0.729	2.5	4.32	400	U H	0	2.5	4.28	U	-		0.4
Carbon 14	15.8	2.2	3.02	50.0	C	16.8	2.3	3.02		6	36	0.5
Nickel 63	<u>-3.91</u>	2.6	4.60	30.0	U NI_L	-1.50	1.9	3.27	U	-		1.5
Technetium 99	4.08	0.36	0.489	12.0	TC	3.98	0.30	0.288		2	33	0.2
Thorium 228	0.494	0.31	0.341	1.00	TH	0.832	0.36	0.366		51	109	1.4
Thorium 230	0.522	0.31	0.235	1.00	TH	0.473	0.30	0.226		10	131	0.2
Thorium 232	0.645	0.31	0.235	1.00	TH	0.769	0.30	0.226		18	93	0.6
Iodine 129	-0.296	0.68	1.56	2.00	U I	0.090	0.59	1.35	U	-		0.9
Radium 226	0.558	0.098	0.078	0.100	RA	0.608	0.10	0.084		9	50	0.5
Radium 228	1.66	0.44	0.989	1.00	AC	1.65	0.44	0.999		1	61	0

200TW1 OU Charactrzn 299E33342-Soil

QC-DUP#1 65615

DUPLICATES

Page 1

SUMMARY DATA SECTION

Page 10

Lab id <u>EBRLNE</u>
Protocol <u>Fluor</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DUP</u>
Version <u>3.06</u>
Report date <u>06/13/08</u>

00000014

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP H3702

7087-001

B1TNR0

D A T A S H E E T

SDG <u>7087</u>	Client/Case no <u>Hanford</u>	<u>SDG_H3702</u>
Contact <u>Melissa C. Mannion</u>	Contract <u>No. 33677</u>	
Lab sample id <u>R804167-01</u>	Client sample id <u>B1TNR0</u>	
Dept sample id <u>7087-001</u>	Location/Matrix <u>C5857, I-098</u>	<u>SOLID</u>
Received <u>04/29/08</u>	Collected/Weight <u>04/01/08 08:20</u>	<u>157 g</u>
% solids <u>89.3</u>	Custody/SAF No <u>F08-070-028</u>	<u>F08-070</u>

ANALYTE	CAS NO	RESULT pCi/g	2σ ERR (COUNT)	MDA pCi/g	RDL pCi/g	QUALI- FIERS	TEST
Tritium	10028-17-8	0	2.5	4.28	400	U	H
Carbon 14	14762-75-5	16.8	2.3	3.02	50.0		C
Nickel 63	13981-37-8	-1.50	1.9	3.27	30.0	U	NI_L
Technetium 99	14133-76-7	3.98	0.30	0.288	12.0		TC
Thorium 228	14274-82-9	0.832	0.36	0.366	1.00		TH
Thorium 230	14269-63-7	0.473	0.30	0.226	1.00		TH
Thorium 232	TH-232	0.769	0.30	0.226	1.00		TH
Iodine 129	15046-84-1	0.090	0.59	1.35	2.00	U	I
Radium 226	13982-63-3	0.608	0.10	0.084	0.100		RA
Radium 228	15262-20-1	1.65	0.44	0.999	1.00		AC

200TW1 OU Charactrzn 299E33342-Soil

DATA SHEETS

Page 1

SUMMARY DATA SECTION

Page 11

Lab id <u>EBRLNE</u>
Protocol <u>Fluor</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>06/13/08</u>

00000015

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test TH Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

Client Hanford
 Contract No. 33677
 Contract SDG H3702

LAB METHOD SUMMARY

THORIUM, ISOTOPIC IN SOLIDS

ALPHA SPECTROSCOPY

RESULTS

LAB RAW SUF-
 SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Thorium 230

Preparation batch 6150-050

R804167-01	7087-001	BITNR0	0.473
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	ok

Nominal values and limits from method RDLs (pCi/g) 1.00
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB	RAW	SUF-	MAX MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD PREPARED	YZED DETECTOR

Preparation batch 6150-050	2σ prep error 8.0 %	Reference Lab Notebook #6150, pg. 55											
R804167-01	BITNR0	0.366	0.250	90	174	63	06/02/08	06/03	SS-061				
R804167-02	Lab Control Sample	0.214	0.250	87	186		06/02/08	06/05	SS-037				
R804167-03	Method Blank	0.418	0.250	84	186		06/02/08	06/05	SS-038				
R804167-04	Duplicate (R804167-01)	0.341	0.250	86	174	63	06/02/08	06/03	SS-064				

Nominal values and limits from method 1.00 0.250 20-105 150 180

PROCEDURES	REFERENCE	THISO_IE_PLATE_AEA
SPP-061	Determination of Moisture Content in Solid Samples	rev 0
SPP-071	Soil Dissolution, > 1.0g Aliquot, rev 5	
CP-900	Thorium in Water and Dissolved Solid Samples by	Extraction Chromatography, rev 1
CP-008	Heavy Element Electroplating, rev 9	

AVERAGES ± 2 SD	MDA	<u>0.335</u> ± <u>0.173</u>
FOR 4 SAMPLES	YIELD	<u>87</u> ± <u>5</u>

METHOD SUMMARIES

Page 1

SUMMARY DATA SECTION

Page 12

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000016

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test AC Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

RADIUM 228 IN SOLIDS

BETA COUNTING

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB RAW SUP-
 SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Radium 228

Preparation batch 6150-050

R804167-01	7087-001	BITNR0	1.65
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	ok

Nominal values and limits from method RDLs (pCi/g) 1.00
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB RAW SUP- MDA ALIQ PREP DILU- YIELD EFF COUNT FWHM DRIFT DAYS ANAL-
 SAMPLE ID TEST FIX CLIENT SAMPLE ID pCi/g g FAC TION % % min keV KeV HELD PREPARED YZED DETECTOR

Preparation batch 6150-050 2σ prep error 10.4 % Reference Lab Notebook #6150, pg. 55

R804167-01	BITNR0	0.999	1.00	75	150	69	06/09/08	06/09	GRB-221
R804167-02	Lab Control Sample	0.967	1.00	73	150		06/09/08	06/09	GRB-222
R804167-03	Method Blank	0.978	1.00	74	150		06/09/08	06/09	GRB-223
R804167-04	Duplicate (R804167-01)	0.989	1.00	74	150	69	06/09/08	06/09	GRB-224

Nominal values and limits from method 1.00 1.00 100 180

PROCEDURES REFERENCE RAISO_SEP_GPC
 SPP-071 Soil Dissolution, > 1.0g Aliquot, rev 5
 DWP-894 Sequential Separation of Actinium-228 and Radium-226 in Drinking Water (>1 Liter Aliquot), rev 0

AVERAGES ± 2 SD MDA 0.983 ± 0.028
 FOR 4 SAMPLES YIELD 74 ± 2

METHOD SUMMARIES

Page 2

SUMMARY DATA SECTION

Page 13

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000017

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test TC Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

TECHNETIUM 99 IN SOLIDS
 BETA COUNTING

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB	RAW	SUF-		Technetium
SAMPLE ID	TEST FIX	PLANCHET	CLIENT SAMPLE ID	99

Preparation batch 6150-050

R804167-01	7087-001	B1TNR0	3.98
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	ok

Nominal values and limits from method RDLs (pCi/g) 12.0
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-		
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD	PREPARED	YZED	DETECTOR

Preparation batch 6150-050 2σ prep error 13.2 % Reference Lab Notebook #6150, pg. 55

R804167-01	B1TNR0	0.288	1.00		91	178	58	05/22/08	05/29	GRB-230
R804167-02	Lab Control Sample	0.542	1.00		100	50		05/22/08	05/27	GRB-232
R804167-03	Method Blank	0.486	1.00		<u>106</u>	50		05/22/08	05/27	GRB-222
R804167-04	Duplicate (R804167-01)	0.489	1.00		58	178	58	05/22/08	05/29	GRB-231

Nominal values and limits from method 12.0 1.00 20-105 50 180

PROCEDURES	REFERENCE	TC99_TR_SEP_GPC
	SPP-062	Sample Aliquoting, rev 0
	CP-431	Technetium-99 Purification of Soil or Resin by Extraction Chromatography, rev 2
	CP-008	Heavy Element Electroplating, rev 9

AVERAGES ± 2 SD	MDA <u>0.451</u> ± <u>0.224</u>
FOR 4 SAMPLES	YIELD <u>89</u> ± <u>43</u>

METHOD SUMMARIES

Page 3

SUMMARY DATA SECTION

Page 14

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000018

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test I Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

IODINE 129 IN SOLIDS
 GAMMA SPECTROSCOPY

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB	RAW	SUF-		
SAMPLE ID	TEST FIX	PLANCHET	CLIENT SAMPLE ID	Iodine 129
Preparation batch 6150-050				
R804167-01		7087-001	B1TNR0	U
R804167-02		7087-002	Lab Control Sample	ok
R804167-03		7087-003	Method Blank	U
R804167-04		7087-004	Duplicate (R804167-01)	- U
Nominal values and limits from method RDLs (pCi/g) 2.00				
200TW1 OU Charactrzn 299E33342-Soil				

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-		
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD	PREPARED	YZED	DETECTOR
Preparation batch 6150-050 2σ prep error 19.4 % Reference Lab Notebook #6150, pg. 55															
R804167-01		B1TNR0	1.35	1.00			59					65	06/05/08	06/05	XSPEC-004
R804167-02		Lab Control Sample	1.88	1.00			81						06/05/08	06/05	XSPEC-002
R804167-03		Method Blank	1.32	1.00			78						06/05/08	06/06	XSPEC-004
R804167-04		Duplicate (R804167-01)	1.56	1.00			68						06/05/08	06/06	XSPEC-002
Nominal values and limits from method 2.00 1.00 20-105 300 180															

PROCEDURES	REFERENCE	I129_SEP_LEPS_GS
SPP-062	Sample Aliquoting, rev 0	
SPP-060	Soil Preparation, rev 0	
CP-024	Iodine-129, Sample Dissolution, rev 5	
CP-530	Iodine-129 Purification, rev 1	

AVERAGES ± 2 SD	MDA	<u>1.53</u>	±	<u>0.516</u>
FOR 4 SAMPLES	YIELD	<u>72</u>	±	<u>20</u>

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000019

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

LAB METHOD SUMMARY

CARBON 14 IN SOLIDS

LIQUID SCINTILLATION COUNTING

Test C Matrix SOLID
SDG 7087
Contact Melissa C. Mannion

Client Hanford
Contract No. 33677
Contract SDG H3702

RESULTS

LAB RAW SUF-
SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Carbon 14

Preparation batch 6150-050

R804167-01	7087-001	B1TNR0	16.8
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	ok

Nominal values and limits from method RDLs (pCi/g) 50.0
200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD PREPARED	YZED DETECTOR

Preparation batch 6150-050 2σ prep error 10.0 % Reference Lab Notebook #6150, pg. 55

R804167-01		B1TNR0	3.02	0.458			100		50			62	05/29/08 06/02 LSC-004
R804167-02		Lab Control Sample	3.46	0.400			100		50				05/29/08 06/02 LSC-004
R804167-03		Method Blank	3.75	0.400			100		50				05/29/08 06/02 LSC-004
R804167-04		Duplicate (R804167-01)	3.02	0.456			100		50			62	05/29/08 06/02 LSC-004

Nominal values and limits from method 50.0 0.400 10 180

PROCEDURES REFERENCE C14_COX_LSC
CP-251 Tritium/Carbon-14 Oxidation, rev 8

AVERAGES ± 2 SD MDA 3.31 ± 0.716
FOR 4 SAMPLES YIELD 100 ± 0

METHOD SUMMARIES

Page 5

SUMMARY DATA SECTION

Page 16

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-LMS
Version 3.06
Report date 06/13/08

00000020

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test H Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

TRITIUM IN SOLIDS

LIQUID SCINTILLATION COUNTING

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB	RAW	SUF-		
SAMPLE ID	TEST FIX	PLANCHET	CLIENT SAMPLE ID	Tritium
Preparation batch 6150-050				
R804167-01		7087-001	BITNR0	U
R804167-02		7087-002	Lab Control Sample	ok
R804167-03		7087-003	Method Blank	U
R804167-04		7087-004	Duplicate (R804167-01)	- U

Nominal values and limits from method RDLs (pCi/g) 400
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-		
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD	PREPARED	YZED	DETECTOR
Preparation batch 6150-050			2σ prep error 10.0 % Reference Lab Notebook #6150, pg. 55												
R804167-01		BITNR0	4.28	0.458			100		50		59	05/29/08	05/30	LSC-006	
R804167-02		Lab Control Sample	5.25	0.400			100		50			05/29/08	05/30	LSC-006	
R804167-03		Method Blank	4.78	0.400			100		50			05/29/08	05/30	LSC-006	
R804167-04		Duplicate (R804167-01)	4.32	0.456			100		50		59	05/29/08	05/30	LSC-006	

Nominal values and limits from method 400 0.400 25 180

PROCEDURES REFERENCE TRITIUM_COX_LSC
 CP-251 Tritium/Carbon-14 Oxidation, rev 8

AVERAGES ± 2 SD MDA 4.66 ± 0.911
 FOR 4 SAMPLES YIELD 100 ± 0

METHOD SUMMARIES

Page 6

SUMMARY DATA SECTION

Page 17

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000021

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test NI L Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

NICKEL 63 IN SOLIDS

LIQUID SCINTILLATION COUNTING

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB RAW SUF-
 SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Nickel 63

Preparation batch 6150-050

R804167-01	7087-001	BITNRO	U
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	- U

Nominal values and limits from method RDLs (pCi/g) 30.0
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/g	g	FAC	TION	%	%	min	keV	KeV	HELD PREPARED	YZED DETECTOR

Preparation batch 6150-050 2σ prep error 11.2 % Reference Lab Notebook #6150, pg. 55

R804167-01	BITNRO		3.27	0.500			92		50		57	05/28/08	05/28	LSC-004
R804167-02	Lab Control Sample		3.15	0.500			96		50			05/28/08	05/28	LSC-004
R804167-03	Method Blank		3.16	0.500			96		50			05/28/08	05/28	LSC-004
R804167-04	Duplicate (R804167-01)		4.60	0.500			67		50		57	05/28/08	05/28	LSC-004

Nominal values and limits from method 30.0 0.500 30-105 25 180

PROCEDURES	REFERENCE	NI63_LSC
	SPP-071	Soil Dissolution, > 1.0g Aliquot, rev 5
	CP-280	Nickel-63 Purification, rev 3

AVERAGES ± 2 SD	MDA	<u>3.54</u> ± <u>1.41</u>
FOR 4 SAMPLES	YIELD	<u>88</u> ± <u>28</u>

METHOD SUMMARIES

Page 7

SUMMARY DATA SECTION

Page 18

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000022

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H3702

Test RA Matrix SOLID
 SDG 7087
 Contact Melissa C. Mannion

LAB METHOD SUMMARY

RADIUM 226 IN SOLIDS
 RADON COUNTING

Client Hanford
 Contract No. 33677
 Contract SDG H3702

RESULTS

LAB RAW SUF-
 SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Radium 226

Preparation batch 6150-050

R804167-01	7087-001	B1TNR0	0.608
R804167-02	7087-002	Lab Control Sample	ok
R804167-03	7087-003	Method Blank	U
R804167-04	7087-004	Duplicate (R804167-01)	ok

Nominal values and limits from method RDLs (pCi/g) 0.100
 200TW1 OU Charactrzn 299E33342-Soil

METHOD PERFORMANCE

LAB RAW SUF- MDA ALIQ PREP DILU- YIELD EFF COUNT FWHM DRIFT DAYS ANAL-
 SAMPLE ID TEST FIX CLIENT SAMPLE ID pCi/g g FAC TION % % min keV KeV HELD PREPARED YZED DETECTOR

Preparation batch 6150-050 2σ prep error 16.4 % Reference Lab Notebook #6150, pg. 55

R804167-01	B1TNR0	0.084	1.00	100	73	56	05/27/08	05/27	RN-011
R804167-02	Lab Control Sample	0.082	1.00	100	73		05/27/08	05/27	RN-012
R804167-03	Method Blank	0.071	1.00	100	73		05/27/08	05/27	RN-013
R804167-04	Duplicate (R804167-01)	0.078	1.00	100	73	56	05/27/08	05/27	RN-014

Nominal values and limits from method 0.100 1.00 20-105 50 180

PROCEDURES REFERENCE 903.1_RA226_LUC
 SPP-071 Soil Dissolution, > 1.0g Aliquot, rev 5
 CP-881 Radium-226 in Water and Dissolved Samples, rev 1

AVERAGES ± 2 SD MDA 0.079 ± 0.011
 FOR 4 SAMPLES YIELD 100 ± 0

METHOD SUMMARIES

Page 8

SUMMARY DATA SECTION

Page 19

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-LMS
 Version 3.06
 Report date 06/13/08

00000023

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087

Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford

Contract No. 33677

Case no SDG_H3702

SAMPLE SUMMARY

The Sample and QC Summary Reports show all samples, including QC samples, reported in one Sample Delivery Group (SDG).

The Sample Summary Report fully identifies client samples and gives the corresponding lab sample identification. The QC Summary Report shows at the sample level how the lab organized the samples into batches and generated QC samples. The Preparation Batch and Method Summary Reports show this at the analysis level.

The following notes apply to these reports:

- * LAB SAMPLE ID is the lab's primary identification for a sample.
- * DEPARTMENT SAMPLE ID is an alternate lab id, for example one assigned by a radiochemistry department in a lab.
- * CLIENT SAMPLE ID is the client's primary identification for a sample. It includes any sample preparation done by the client that is necessary to identify the sample.
- * QC BATCH is a lab assigned code that groups samples to be processed and QCed together. These samples should have similar matrices.

QC BATCH is not necessarily the same as SDG, which reflects samples received and reported together.

- * All Lab Control Samples, Method Blanks, Duplicates and Matrix Spikes are shown that QC any of the samples. Due to possible reanalyses, not all results for all these QC samples may be relevant to the SDG. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 20

Lab id EBRLNE

Protocol Fluor

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 06/13/08

00000024

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG H3702

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches in one Sample Delivery Group (SDG) with information necessary to check the completeness and consistency of the SDG.

The following notes apply to this report:

- * The preparation batches are shown in the same order as the Method Summary Reports are printed.
- * Only analyses of planchets relevant to the SDG are included.
- * Each preparation batch should have at least one Method Blank and LCS in it to validate client sample results.
- * The QUALIFIERS shown are all qualifiers other than U, J, B, L and H that occur on any analysis in the preparation batch. The Method Summary Report has these qualifiers on a per sample basis.

These qualifiers should be reviewed as follows:

- X Some data has been manually entered or modified. Transcription errors are possible.
- P One or more results are 'preliminary'. The data is not ready for final reporting.
- 2 There were two or more results for one analyte on one planchet imported at one time. The results in DVD may not be the same as on the raw data sheets.

Other lab defined qualifiers may occur. In general, these should be addressed in the SDG narrative.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 21

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000025

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG_H3702

WORK SUMMARY

The Work Summary Report shows all samples, including QC samples, and all relevant analyses in one Sample Delivery Group (SDG). This report is often useful as supporting documentation for an invoice.

The following notes apply to this report:

- * TEST is a code for the method used to measure associated analytes. Results and related information for each analyte are on the Data Sheet Report. In special cases, a test code used in the summary data section is not the same as in associated raw data. In this case, both codes are shown on the Work Summary.
- * SUFFIX is the lab's code to distinguish multiple analyses (recounts, reworks, reanalyses) of a fraction of the sample. The suffix indicates which result is being reported. An empty suffix normally identifies the first attempt to analyze the sample.
- * The LAB SAMPLE ID, TEST and SUFFIX uniquely identify all supporting data for a result. The Method Summary Report for each TEST has method performance data, such as yield, for each lab sample id and suffix and procedures used in the method.
- * PLANCHET is an alternate lab identifier for work done for one test. It, combined with the TEST and SUFFIX, may be the best link to raw data.
- * For QC samples, only analyses that directly QC some regular sample are shown. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.
- * The SAS (Special Analytical Services) Number is a client or lab assigned code that reflects special processing for samples, such as rapid turn around. Counts of tests done are lists by SAS number since it is likely to affect prices.

REPORT GUIDES

Page 3

SUMMARY DATA SECTION

Page 22

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000026

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG H3702

DATA SHEET

The Data Sheet Report shows all results and primary supporting information for one client sample or Method Blank. This report corresponds to both the CLP Inorganics and Organics Data Sheet.

The following notes apply to this report:

- * TEST is a code for the method used to measure an analyte. If the TEST is empty, no data is available; the analyte was not analyzed for.
- * The LAB SAMPLE ID and TEST uniquely identify work within the Summary Data Section of a Data Package. The Work Summary and Method Summary Reports further identify raw data that underlies this work.

The Method Summary Report for each TEST has method performance data, such as yield, for each Lab Sample ID and a list of procedures used in the method.

- * ERRORS can be labeled TOTAL or COUNT. TOTAL implies a preparation (non-counting method) error has been added, as square root of sum of squares, to the counting error denoted by COUNT. The preparation errors, which may vary by preparation batch, are shown on the Method Summary Report.
- * A RESULT can be 'N.R.' (Not Reported). This means the lab did this work but chooses not to report it now, possibly because it was reported at another time.
- * When reporting a Method Blank, a RESULT can be 'N.A.' (Not Applicable). This means there is no reported client sample work in the same preparation batch as the Blank's result. This is likely to occur when the Method Blank is associated with reanalyses of selected work for a few samples in the SDG.

The following qualifiers are defined by the DVD system:

U The RESULT is less than the MDA (Minimum Detectable Activity).

REPORT GUIDES

Page 4

SUMMARY DATA SECTION

Page 23

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000027

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
 Contract No. 33677
 Case no SDG H3702

DATA SHEET

If the MDA is blank, the ERROR is used as the limit.

- J The RESULT is less than the RDL (Required Detection Limit) and no U qualifier is assigned.
 - B A Method Blank associated with this sample had a result without a U flag and, after correcting for possibly different aliquots, that result is greater than or equal to the MDA for this sample.
- Normally, B is not assigned if U is. When method blank subtraction is shown on this report, B flags are assigned based on the unsubtracted values while U's are assigned based on the subtracted ones. Both flags can be assigned in this case.
- For each sample result, all Method Blank results in the same preparation batch are compared. The Method Summary Report documents this and other QC relationships.
- L Some Lab Control Sample that QC's this sample had a low recovery. The lab can disable assignment of this qualifier.
 - H Similar to 'L' except the recovery was high.
 - P The RESULT is 'preliminary'.
 - X Some data necessary to compute the RESULT, ERROR or MDA was manually entered or modified.
 - 2 There were two or more results available for this analyte. The reported result may not be the same as in the raw data.

Other qualifiers are lab defined. Definitions should be in the SDG narrative.

The following values are underlined to indicate possible problems:

- * An MDA is underlined if it is bigger than its RDL.

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 06/13/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087

Contact Melissa C. Mannion

Client Hanford

Contract No. 33677

Case no SDG_H3702

GUIDE , cont .

DATA SHEET

- * An ERROR is underlined if the 1.645 sigma counting error is bigger than both the MDA and the RESULT, implying that the MDA may not be a good estimate of the 'real' minimum detectable activity.
- * A negative RESULT is underlined if it is less than the negative of its 2 sigma counting ERROR.
- * When reporting a Method Blank, a RESULT is underlined if greater than its MDA. If the MDA is blank, the 2 sigma counting error is used in the comparison.

REPORT GUIDES

Page 6

SUMMARY DATA SECTION

Page 25

Lab id EBRLNE

Protocol Fluor

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 06/13/08

00000029

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
 Contract No. 33677
 Case no SDG H3702

LAB CONTROL SAMPLE

The Lab Control Sample Report shows all results, recoveries and primary supporting information for one Lab Control Sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. Refer to its Report Guide for details.
- * An amount ADDED is the lab's value for the actual amount spiked into this sample with its ERROR an estimate of the error of this amount.

An amount added is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is RESULT divided by ADDED expressed as a percent.
- * The first, computed limits for the recovery reflect:
 1. The error of RESULT, including that introduced by rounding the result prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.
 2. The error of ADDED.
 3. A lab specified, per analyte bias. The bias changes the center of the computed limits.
- * The second limits are protocol defined upper and lower QC limits for the recovery.
- * The recovery is underlined if it is outside either of these ranges.

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 06/13/08

00000030

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG H3702

DUPLICATE

The Duplicate Report shows all results, differences and primary supporting information for one Duplicate and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Duplicate and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Duplicate has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * The RPD (Relative Percent Difference) is the absolute value of the difference of the RESULTS divided by their average expressed as a percent.

If both RESULTS are less than their MDAs, no RPD is computed and a '-' is printed.

For an analyte, if the lab did work for both samples but has data for only one, the MDA from the sample with data is used as the other's result in the RPD.

- * The first, computed limit is the sum, as square root of sum of squares, of the errors of the results divided by the average result as a percent, hence the relative error of the difference rather than the error of the relative difference. The errors include those introduced by rounding the RESULTS prior to printing.

If this limit is labeled TOT, it includes the preparation error in the RESULTS. If labeled CNT, it does not.

This value reported for this limit is at most 999.

- * The second limit for the RPD is the larger of:

1. A fixed percentage specified in the protocol.

REPORT GUIDES

Page 8

SUMMARY DATA SECTION

Page 27

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000031

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. 33677
Case no SDG H3702

DUPLICATE

2. A protocol factor (typically 2) times the average MDA as a percent of the average result. This limit applies when the results are close to the MDAs.

- * The RPD is underlined if it is greater than either limit.
- * If specified by the lab, the second limit column is replaced by the Difference Error Ratio (DER), which is the absolute value of the difference of the results divided by the quadratic sum of their one sigma errors, the same errors as used in the first limit.

Except for differences due to rounding, the DER is the same as the RPD divided by the first RPD limit with the limit scaled to 1 sigma.

- * The DER is underlined if it is greater than the sigma factor, typically 2 or 3, shown in the header for the first RPD limit.

REPORT GUIDES

Page 9

SUMMARY DATA SECTION

Page 28

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000032

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. 33677
Case no SDG H3702

MATRIX SPIKE

The Matrix Spike Report shows all results, recoveries and primary supporting information for one Matrix Spike and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Spiked and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Spike has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * An amount ADDED is the lab's value for the actual amount spiked into the Spike sample with its ERROR an estimate of the error of this amount.

An amount is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is the Spike RESULT minus the Original RESULT divided by ADDED expressed as a percent.

- * The first, computed limits for the recovery reflect:

1. The errors of the two RESULTS, including those introduced by rounding them prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.

2. The error of ADDED.

3. A lab specified, per analyte bias. The bias changes the center of the computed limits.

- * The second limits are protocol defined upper and lower QC limits

REPORT GUIDES

Page 10

SUMMARY DATA SECTION

Page 29

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000033

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087

Contact Melissa C. Mannion

Client Hanford

Contract No. 33677

Case no SDG H3702

GUIDE , c o n t .

MATRIX SPIKE

for the recovery.

These limits are left blank if the Original RESULT is more than a protocol defined factor (typically 4) times ADDED. This is a way of accounting for that when the spike is small compared to the amount in the original sample, the recovery is unreliable.

- * The recovery is underlined (out of spec) if it is outside either of these ranges.

REPORT GUIDES

Page 11

SUMMARY DATA SECTION

Page 30

Lab id EBRLNE

Protocol Fluor

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 06/13/08

00000034

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
 Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
 Contract No. 33677
 Case no SDG H3702

METHOD SUMMARY

The Method Summary Report has two tables. One shows up to five results measured using one method. The other has performance data for the method. There is one report for each TEST, as used on the Data Sheet Report.

The following notes apply to this report:

- * Each table is subdivided into sections, one for each preparation batch. A preparation batch is a group of aliquots prepared at roughly the same time in one work area of the lab using the same method.

There should be Lab Control Sample and Method Blank results in each preparation batch since this close correspondence makes the QC meaningful. Depending on lab policy, Duplicates need not occur in each batch since they QC sample dependencies such as matrix effects.

- * The RAW TEST column shows the test code used in the raw data to identify a particular analysis if it is different than the test code in the header of the report. This occurs in special cases due to method specific details about how the lab labels work.

The Lab Sample or Planchet ID combined with the (Raw) Test Code and Suffix uniquely identify the raw data for each analysis.

- * If a result is less than both its MDA and RDL, it is replaced by just 'U' on this report. If it is greater than or equal to the RDL but less than the MDA, the result is shown with a 'U' flag.

The J and X flags are as on the data sheet.

- * Non-U results for Method Blanks are underlined to indicate possible contamination of other samples in the preparation batch. The Method Blank Report has supporting data.

- * Lab Control Sample and Matrix Spike results are shown as: ok, No data, LOW or HIGH, with the last two underlined. 'No data'

Lab id EBRLNE
 Protocol Fluor
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 06/13/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

GUIDE , c o n t .

Client Hanford
Contract No. 33677
Case no SDG_H3702

METHOD SUMMARY

means no amount ADDED was specified. 'LOW' and 'HIGH' correspond to when the recovery is underlined on the Lab Control Sample or Matrix Spike Report. See these reports for supporting data.

- * Duplicate sample results are shown as: ok, No data, or OUT, with the last two underlined. 'No data' means there was no original sample data found for this duplicate. 'OUT' corresponds to when the RPD is underlined on the Duplicate Report. See this report for supporting data.
- * If the MDA column is labeled 'MAX MDA', there was more than one result measured by the reported method and the MDA shown is the largest MDA. If not all these results have the same RDL, the MAX MDA reflects only those results with RDL equal to the smallest one.

MDAs are underlined if greater than the printed RDL.
- * Aliquots are underlined if less than the nominal value specified for the method.
- * Preparation factors are underlined if greater than the nominal value specified for the method.
- * Dilution factors are underlined if greater than the nominal value specified for the method.
- * Residues are underlined if outside the range specified for the method. Residues are not printed if yields are.
- * Yields, which may be gravimetric, radiometric or some type of recovery depending on the method, are underlined if outside the range specified for the method.
- * Efficiencies are underlined if outside the range specified for the method. Efficiencies are detector and geometry dependent so this test is only approximate.

REPORT GUIDES

Page 13

SUMMARY DATA SECTION

Page 32

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000036

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. 33677
Case no SDG H3702

METHOD SUMMARY

- * Count times are underlined if less than the nominal value specified for the method.
- * Resolutions (as FWHM; Full Width at Half Max) are underlined if greater than the method specified limit.
- * Tracer drifts are underlined if their absolute values are greater than the method specified limit. Tracer drifts are not printed if percent moistures are.
- * Days Held are underlined if greater than the holding time specified in the protocol.
- * Analysis dates are underlined if before their planchet's preparation date or, if a limit is specified, too far after it.

For some methods, ratios as percentages and error estimates for them are computed for pairs of results. A ratio column header like '1+3' means the ratio of the first result column and the third result column.

Ratios are not computed for Lab Control Sample, Method Blank or Matrix Spike results since their matrices are not necessarily similar to client samples'.

The error estimate for a ratio of results from one planchet reflects only counting errors since other errors should be correlated. For a ratio involving different planchets, if QC limits are computed based on total errors, the error for the ratio allows for the preparation errors for the planchets.

The ratio is underlined (out of spec) if the absolute value of its difference from the nominal value is greater than its error estimate. If no nominal value is specified, this test is not done.

For Gross Alpha or Gross Beta results, there may be a column showing the sum of other Alpha or Beta emitters. This sum includes all relevant

REPORT GUIDES

Page 14

SUMMARY DATA SECTION

Page 33

Lab id EBRLNE
Protocol Fluor
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 06/13/08

00000037

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H3702

SDG 7087

Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford

Contract No. 33677

Case no SDG H3702

METHOD SUMMARY

results in the DVD database, whether reported or not. Results in the sum are weighted by a particles/decay value specified by the lab for each relevant analyte. Results less than their MDA are not included. No sums are computed for Lab Control, Method Blank or Matrix Spike samples since their various planchets may not be physically related.

If a ratio of total isotopic to Gross Alpha or Beta is shown, the error for the ratio reflects both the error in the Gross result and the sum, as square root of sum of squares, of the errors in the isotopic results.

For total elemental uranium or thorium results, there may be a column showing the total weight computed from associated isotopic results. Ignoring results less than their MDAs, this is a weighted sum of the isotopic results. The weights depend on the molecular weight and half-life of each isotope so as to convert activities (decays) to weight (atoms).

If a ratio of total computed to measured elemental uranium or thorium is shown, the error for the ratio reflects the errors in all the measurements.

REPORT GUIDES

Page 15

SUMMARY DATA SECTION

Page 34

Lab id EBRLNE

Protocol Fluor

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 06/13/08

00000038

COLLECTOR

NCO SAMPLER

COMPANY CONTACT

TRENT, STEVE

TELEPHONE NO.

373-5869

PROJECT COORDINATOR

WIDRIG, DL

PRICE CODE 8N

DATA TURNAROUND

SAMPLING LOCATION

C5857, I-098

2/c 4-16-08

PROJECT DESIGNATION

200-TW-1 OU Characterization for Well 299-E33-342 - Soil

H3702 (7087)

SAF NO.

F08-070

AIR QUALITY

45 Days / 45 Days

ICE CHEST NO.

GLSS-6

FIELD LOGBOOK NO.

HNF-N-4881 Pg 12

ACTUAL SAMPLE DEPTH

1898-192.3

COA

123163ES20

METHOD OF SHIPMENT

FEDERAL EXPRESS

SHIPPED TO

Eberline Services

OFFSITE PROPERTY NO.

See PTR

BILL OF LADING/AIR BILL NO.

See PTR *7984 2148 6655*

MATRIX*

- A=Air
- DL=Drum
- Liquids
- DS=Drum
- Solids
- L=Liquid
- O=Oil
- S=Soil
- SE=Sediment
- T=Tissue
- V=Vegetation
- W=Water
- WI=Wipe
- X=Other

POSSIBLE SAMPLE HAZARDS/ REMARKS

Contains Radioactive Material at concentrations that are not regulated for transportation per 49 CFR but are not releasable per DOE Order 5400.5 (1990/1993)

PRESERVATION

None

TYPE OF CONTAINER

G/P

NO. OF CONTAINER(S)

1

VOLUME

120mL

SAMPLE ANALYSIS

SEE ITEM (1) IN SPECIAL INSTRUCTIONS

SPECIAL HANDLING AND/OR STORAGE

Radioactive Tie To B1TNP7

SAMPLE NO.	MATRIX*	SAMPLE DATE	SAMPLE TIME
B1TNR0	SOIL	APR 01 2008	0820 ✓

024875

CHAIN OF POSSESSION

SIGN/ PRINT NAMES

SPECIAL INSTRUCTIONS

RELINQUISHED BY/REMOVED FROM *[Signature]* DATE/TIME *APR 01 2008*

RECEIVED BY/STORED IN *MO745 F08-070 #1* DATE/TIME *APR 01 2008 1550*

** The 200 Area S&GRP Characterization and Monitoring Sampling and Analysis GKI applies to this SAF. (1)Carbon-14; Technetium-99 {Technetium-99} Iodine-129; TRITIUM - MIDLEVEL; Isotopic Thorium; Radium -226 {Radium-226} Nickel-63; Radium -228 {Radium-228}

RELINQUISHED BY/REMOVED FROM *170 745 REF #1* DATE/TIME *4-16-08 0700*

RECEIVED BY/STORED IN *El Kawan / E. J. ... 4-16-08 0700* DATE/TIME

RELINQUISHED BY/REMOVED FROM *Ed Kawan* DATE/TIME *4-16-08 0700*

RECEIVED BY/STORED IN *Fed X* DATE/TIME

RELINQUISHED BY/REMOVED FROM *Fed Ex* DATE/TIME *4-17-08 0935*

RECEIVED BY/STORED IN *[Signature]* DATE/TIME *4-17-08 0935*

RELINQUISHED BY/REMOVED FROM *[Signature]* DATE/TIME *4-22-08 1120*

RECEIVED BY/STORED IN *[Signature]* DATE/TIME

RELINQUISHED BY/REMOVED FROM *UPS* DATE/TIME

RECEIVED BY/STORED IN *[Signature]* DATE/TIME *04/24/08 10:30*

RELINQUISHED BY/REMOVED FROM

RECEIVED BY/STORED IN

LABORATORY SECTION

RECEIVED BY

TITLE

DATE/TIME

FINAL SAMPLE DISPOSITION

DISPOSAL METHOD

DISPOSED BY

DATE/TIME

000000039



RICHMOND, CA LABORATORY

SAMPLE RECEIPT CHECKLIST

4-29/08

Client: F. HANFORD / LIONVILLE City TUCKER ANDERSON / OXTON State WA / PA

Date/Time received 04/29/08 10:30 CoC No. FC8070-028

Container I.D. No. Guss - 6 (Box) Requested TAT (Days) 45 P.O. Received Yes [] No []

INSPECTION

- 1. Custody seals on shipping container intact? Yes [] No [] N/A []
 - 2. Custody seals on shipping container dated & signed? Yes [] No [] N/A []
 - 3. Custody seals on sample containers intact? Yes [] No [] N/A []
 - 4. Custody seals on sample containers dated & signed? Yes [] No [] N/A []
 - 5. Packing material is: Wet [] Dry []
 - 6. Number of samples in shipping container 1 Sample Matrix S
 - 7. Number of containers per sample: 1 (Or see CoC)
 - 8. Samples are in correct container Yes [] No []
 - 9. Paperwork agrees with samples? Yes [] No []
 - 10. Samples have Tape [] Hazard labels [] Rad labels [] Appropriate sample labels []
 - 11. Samples are In good condition [] Leaking [] Broken Container [] Missing []
 - 12. Samples are Preserved [] Not preserved [] pH Preservative
 - 13. Describe any anomalies:
14. Was P.M. notified of any anomalies? Yes [] No [] Date
15. Inspected by AFK Date 04/29/08 Time 10:45

Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	Wipe	Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	wipe
<u>BITNRG</u>	<u><60</u>						

Ion Chamber Ser. No. Calibration date

Alpha Meter Ser. No. Calibration date

Beta/Gamma Meter Ser. No. 100482 Calibration date 09/24/07

00000040