

RECEIVED OCTOBER 04, 2010

0094492



**EBERLINE**  
SERVICES

EBERLINE ANALYTICAL CORPORATION  
2030 Wright Avenue  
Richmond, California 94804-3845  
Phone (510) 235-2633 Fax (510) 235-0438  
Toll Free (800) 841-5487  
www.eberlineservices.com

October 4, 2010

Mr. Michael Neely  
CH2M Hill Plateau Remediation Company  
P.O. Box 1600  
Mail Stop – B6-06  
Richland, WA 99352

Reference: **P.O. #33677**  
**Eberline Analytical S0-09-128-7899, SDG H4384**

Dear Mr. Neely:

Enclosed is a data report for two water samples designated under SAF No. F10-275 received at Eberline Analytical on September 14, 2010. The samples were analyzed according to the accompanying chain-of-custody documents.


Please call if you have any questions concerning this report.

Sincerely,

*for MWT*  
N. Joseph Verville  
Client Services Manager

MWT/ljb

Enclosure: Data Package

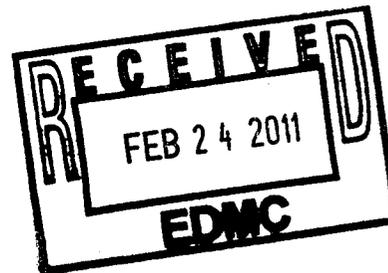



EBRLINE SERVICES / RICHMOND  
SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

S U M M A R Y   D A T A   S E C T I O N

T A B L E   O F   C O N T E N T S				
About this section	.	.	.	1
Sample Summaries	.	.	.	3
Prep Batch Summary	.	.	.	5
Work Summary	.	.	.	6
Method Blanks	.	.	.	7
Lab Control Samples	.	.	.	8
Duplicates	.	.	.	9
Data Sheets	.	.	.	10
Method Summaries	.	.	.	12
Report Guides	.	.	.	13
End of Section	.	.	.	27

VB

Prepared by



Reviewed by

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-TOC  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

ABOUT THE DATA SUMMARY SECTION

The Data Summary Section of a Data Package has all data, in several useful orders, necessary for first level, routine review of the data package for a Sample Delivery Group (SDG). This section follows the Data Package Narrative, which has an overview of the data package and a discussion of special problems. It is followed by the Raw Data Section, which has full details.

The Data Summary Section has several groups of reports:

**SAMPLE SUMMARIES**

The Sample and QC Summary Reports show all samples, including QC samples, reported in one SDG. These reports cross-reference client and lab sample identifiers.

**PREPARATION BATCH SUMMARY**

The Preparation Batch Summary Report shows all preparation batches (lab groupings reflecting how work was organized) relevant to the reported SDG with information necessary to check the completeness and consistency of the SDG.

**WORK SUMMARY**

The Work Summary Report shows all samples and work done on them relevant to the reported SDG.

**METHOD BLANKS**

The Method Blank Reports, one for each Method Blank relevant to the SDG, show all results and primary supporting information for the blanks.

**LAB CONTROL SAMPLES**

The Lab Control Sample Reports, one for each Lab Control Sample relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 1

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

ABOUT THE DATA SUMMARY SECTION

DUPLICATES

The Duplicate Reports, one for each Duplicate and Original sample pair relevant to the SDG, show all results, differences and primary supporting information for these QC samples.

MATRIX SPIKES

The Matrix Spike Reports, one for each Spiked and Original sample pair relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

DATA SHEETS

The Data Sheet Reports, one for each client sample in the SDG, show all results and primary supporting information for these samples.

METHOD SUMMARIES

The Method Summary Reports, one for each test used in the SDG, show all results, QC and method performance data for one analyte on one or two pages. (A test is a short code for the method used to do certain work to the client's specification.)

REPORT GUIDES

The Report Guides, one for each of the above groups of reports, have documentation on how to read the associated reports.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 2

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

**EBERLINE SERVICES/RICHMOND**

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

**LAB SAMPLE SUMMARY**

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

LAB SAMPLE ID	CLIENT SAMPLE ID	LOCATION	MATRIX	LEVEL	SAF NO	CHAIN OF CUSTODY	COLLECTED
S009128-01	B28157	216-N-6 EB 1 Day 1	WATER		F10-275	F10-275-005	09/10/10 07:55
S009128-02	B28159	216-N-6 FB	WATER		F10-275	F10-275-007	09/10/10 09:10
S009128-03	Lab Control Sample		WATER		F10-275		
S009128-04	Method Blank		WATER		F10-275		
S009128-05	Duplicate (S009128-01)	216-N-6 EB 1 Day 1	WATER		F10-275		09/10/10 07:55

LAB SUMMARY

Page 1

SUMMARY DATA SECTION

Page 3

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-LS  
 Version 3.06  
 Report date 10/04/10

**EBERLINE SERVICES/RICHMOND**

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

**QC SUMMARY**

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

QC BATCH	CHAIN OF CUSTODY	CLIENT SAMPLE ID	MATRIX	SOLIDS	SAMPLE AMOUNT	BASIS AMOUNT	DAYS SINCE RECEIVED	LAB COLL	DEPARTMENT SAMPLE ID
7899	F10-275-005	B28157	WATER		1.0 L		09/14/10 4	S009128-01	7899-001
	F10-275-007	B28159	WATER		1.0 L		09/14/10 4	S009128-02	7899-002
		Method Blank	WATER					S009128-04	7899-004
		Lab Control Sample	WATER					S009128-03	7899-003
		Duplicate (S009128-01)	WATER		1.0 L		09/14/10 4	S009128-05	7899-005

QC SUMMARY

Page 1

SUMMARY DATA SECTION

Page 4

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-QS  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

PREP BATCH SUMMARY

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

TEST MATRIX	METHOD	PREPARATION ERROR		PLANCHETS ANALYZED					QUALI-	
		BATCH	2σ %	CLIENT	MORE	RE	BLANK	LCS		DUP/ORIG
Liquid Scintillation Counting										
NI_L	WATER	Nickel-63 in Liquid	7270-040	11.2	2		1	1	1/1	

Duplicates and Matrix Spikes are those with original (Client) sample in this Sample Delivery Group.  
 Blank and LCS planchets are those in the same preparation batch as some Client, Duplicate or Spike sample.

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-PBS  
 Version 3.06  
 Report date 10/04/10

**EBERLINE SERVICES/RICHMOND**

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

**LAB WORK SUMMARY**

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

LAB SAMPLE	CLIENT SAMPLE ID				SUF-					
COLLECTED	LOCATION	MATRIX			FIX	ANALYZED	REVIEWED	BY	METHOD	
RECEIVED	CUSTODY	SAF No	PLANCHET	TEST						
S009128-01	B28157		7899-001	NI_L		09/25/10	10/01/10	BW	Nickel-63 in Liquid	
09/10/10	216-N-6 EB 1 Day 1	WATER								
09/14/10	F10-275-005	F10-275								
S009128-02	B28159		7899-002	NI_L		09/25/10	10/01/10	BW	Nickel-63 in Liquid	
09/10/10	216-N-6 FB	WATER								
09/14/10	F10-275-007	F10-275								
S009128-03	Lab Control Sample		7899-003	NI_L		09/25/10	10/01/10	BW	Nickel-63 in Liquid	
		WATER								
		F10-275								
S009128-04	Method Blank		7899-004	NI_L		09/25/10	10/01/10	BW	Nickel-63 in Liquid	
		WATER								
		F10-275								
S009128-05	Duplicate (S009128-01)		7899-005	NI_L		09/25/10	10/01/10	BW	Nickel-63 in Liquid	
09/10/10	216-N-6 EB 1 Day 1	WATER								
09/14/10		F10-275								

**COUNTS OF TESTS BY SAMPLE TYPE**

TEST	SAF No	METHOD	REFERENCE	CLIENT	MORE	RE	BLANK	LCS	DUP	SPIKE	TOTAL
NI_L	F10-275	Nickel-63 in Liquid	NI63_LSC	2			1	1	1		5
TOTALS				2			1	1	1		5

WORK SUMMARY

Page 1

SUMMARY DATA SECTION

Page 6

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-LWS  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

7899-004

Method Blank

METHOD BLANK

SDG <u>7899</u>	Client/Case no <u>CHPRC</u>	<u>SDG H4384</u>
Contact <u>N. Joseph Verville</u>	Contract <u>No. 33677</u>	
Lab sample id <u>S009128-04</u>	Client sample id <u>Method Blank</u>	
Dept sample id <u>7899-004</u>	Material/Matrix <u>WATER</u>	
	SAP No <u>F10-275</u>	

ANALYTE	CAS NO	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST
Nickel 63	13981-37-8	-1.44	1.5	2.65	15.0	U	NI_L

QC-BLANK #75001

METHOD BLANKS

Page 1

SUMMARY DATA SECTION

Page 7

Lab id <u>EBRLNE</u>
Protocol <u>CHPRC</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>10/04/10</u>



**EBERLINE SERVICES/RICHMOND**

SAMPLE DELIVERY GROUP H4384

7899-005

B28157

**DUPLICATE**

SDG <u>7899</u> Contact <u>N. Joseph Verville</u> DUPLICATE Lab sample id <u>S009128-05</u> Dept sample id <u>7899-005</u>	ORIGINAL Lab sample id <u>S009128-01</u> Dept sample id <u>7899-001</u> Received <u>09/14/10</u>	Client/Case no <u>CHPRC</u> <u>SDG H4384</u> Contract No. <u>33677</u> Client sample id <u>B28157</u> Location/Matrix <u>216-N-6 EB 1 Day 1</u> <u>WATER</u> Collected/Volume <u>09/10/10 07:55</u> <u>1.0 L</u> Custody/SAF No <u>F10-275-005</u> <u>F10-275</u>
--	---	--

ANALYTE	DUPLICATE pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST	ORIGINAL pCi/L	2σ ERR (COUNT)	MDA pCi/L	QUALI- FIERS	RPD %	3σ TOT	DER σ
Nickel 63	0.963	1.4	2.38	15.0	U	NI_L	-1.45	1.5	2.67	U	-	2.3	

QC-DUP#1 75002

ARRA 216-N-6 Sampling and Analysis - Verification - QC Sampling

Lab id <u>EBRLNE</u>
Protocol <u>CHPRC</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DUP</u>
Version <u>3.06</u>
Report date <u>10/04/10</u>

**EBERLINE SERVICES / RICHMOND**  
**SAMPLE DELIVERY GROUP H4384**

7899-001

B28157

**DATA SHEET**

SDG <u>7899</u>	Client/Case no <u>CHPRC</u>	SDG <u>H4384</u>
Contact <u>N. Joseph Verville</u>	Contract No. <u>33677</u>	
Lab sample id <u>S009128-01</u>	Client sample id <u>B28157</u>	
Dept sample id <u>7899-001</u>	Location/Matrix <u>216-N-6 EB 1 Day 1</u>	<u>WATER</u>
Received <u>09/14/10</u>	Collected/Volume <u>09/10/10 07:55</u>	<u>1.0 L</u>
	Custody/SAF No <u>F10-275-005</u>	<u>F10-275</u>

ANALYTE	CAS NO	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST
Nickel 63	13981-37-8	-1.45	1.5	2.67	15.0	U	NI_L

ARRA 216-N-6 Sampling and Analysis - Verification - QC  
 Sampling

DATA SHEETS

Page 1

SUMMARY DATA SECTION

Page 10

Lab id <u>EBRLNE</u>
Protocol <u>CHPRC</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>10/04/10</u>

**EBERLINE SERVICES / RICHMOND**  
**SAMPLE DELIVERY GROUP H4384**

7899-002

B28159

**DATA SHEET**

SDG <u>7899</u>	Client/Case no <u>CHPRC</u>	SDG <u>H4384</u>
Contact <u>N. Joseph Verville</u>	Contract <u>No. 33677</u>	
Lab sample id <u>S009128-02</u>	Client sample id <u>B28159</u>	
Dept sample id <u>7899-002</u>	Location/Matrix <u>216-N-6 FB</u>	<u>WATER</u>
Received <u>09/14/10</u>	Collected/Volume <u>09/10/10 09:10</u>	<u>1.0 L</u>
	Custody/SAF No <u>F10-275-007</u>	<u>F10-275</u>

ANALYTE	CAS NO	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST
Nickel 63	13981-37-8	-0.781	1.4	2.54	15.0	U	NI_L

ARRA 216-N-6 Sampling and Analysis - Verification - QC  
 Sampling

Lab id <u>EBRLNE</u>
Protocol <u>CHPRC</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>10/04/10</u>

**EBERLINE SERVICES/RICHMOND**

SAMPLE DELIVERY GROUP H4384

Test NI L Matrix WATER  
 SDG 7899  
 Contact N. Joseph Verville

**LAB METHOD SUMMARY**

NICKEL-63 IN LIQUID

LIQUID SCINTILLATION COUNTING

Client CHPRC  
 Contract No. 33677  
 Contract SDG H4384

**RESULTS**

LAB RAW SUF-  
 SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Nickel 63

Preparation batch 7270-040

S009128-01	7899-001	B28157	U
S009128-02	7899-002	B28159	U
S009128-03	7899-003	Lab Control Sample	ok
S009128-04	7899-004	Method Blank	U
S009128-05	7899-005	Duplicate (S009128-01)	- U

Nominal values and limits from method RDLs (pCi/L) 15.0

**METHOD PERFORMANCE**

LAB RAW SUF- MDA ALIQ PREP DILU- YIELD EFF COUNT FWHM DRIFT DAYS ANAL-  
 SAMPLE ID TEST FIX CLIENT SAMPLE ID pCi/L L FAC TION % % min keV KeV HELD PREPARED YZED DETECTOR

Preparation batch 7270-040 2σ prep error 11.2 % Reference Lab Notebook No. 7270 pg.040

S009128-01	B28157	2.67	0.500	96	50	15	09/24/10	09/25	LSC-005
S009128-02	B28159	2.54	0.500	99	50	15	09/24/10	09/25	LSC-005
S009128-03	Lab Control Sample	2.61	0.500	97	50		09/24/10	09/25	LSC-005
S009128-04	Method Blank	2.65	0.500	96	50		09/24/10	09/25	LSC-005
S009128-05	Duplicate (S009128-01)	2.38	0.500	105	50	15	09/24/10	09/25	LSC-005

Nominal values and limits from method 15.0 0.500 40-110 50 180

PROCEDURES REFERENCE NI63\_LSC  
 SPP-040 Environmental Water Dissolution, rev 2  
 CP-280 Nickel-63 Purification, rev 5

AVERAGES ± 2 SD MDA 2.57 ± 0.235  
 FOR 5 SAMPLES YIELD 99 ± 8

METHOD SUMMARIES

Page 1

SUMMARY DATA SECTION

Page 12

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-LMS  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

SAMPLE SUMMARY

The Sample and QC Summary Reports show all samples, including QC samples, reported in one Sample Delivery Group (SDG).

The Sample Summary Report fully identifies client samples and gives the corresponding lab sample identification. The QC Summary Report shows at the sample level how the lab organized the samples into batches and generated QC samples. The Preparation Batch and Method Summary Reports show this at the analysis level.

The following notes apply to these reports:

- \* LAB SAMPLE ID is the lab's primary identification for a sample.
- \* DEPARTMENT SAMPLE ID is an alternate lab id, for example one assigned by a radiochemistry department in a lab.
- \* CLIENT SAMPLE ID is the client's primary identification for a sample. It includes any sample preparation done by the client that is necessary to identify the sample.
- \* QC BATCH is a lab assigned code that groups samples to be processed and QCed together. These samples should have similar matrices.

QC BATCH is not necessarily the same as SDG, which reflects samples received and reported together.

- \* All Lab Control Samples, Method Blanks, Duplicates and Matrix Spikes are shown that QC any of the samples. Due to possible reanalyses, not all results for all these QC samples may be relevant to the SDG. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 13

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899

Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC

Contract No. 33677

Case no SDG H4384

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches in one Sample Delivery Group (SDG) with information necessary to check the completeness and consistency of the SDG.

The following notes apply to this report:

- \* The preparation batches are shown in the same order as the Method Summary Reports are printed.
- \* Only analyses of planchets relevant to the SDG are included.
- \* Each preparation batch should have at least one Method Blank and LCS in it to validate client sample results.
- \* The QUALIFIERS shown are all qualifiers other than U, J, B, L and H that occur on any analysis in the preparation batch. The Method Summary Report has these qualifiers on a per sample basis.

These qualifiers should be reviewed as follows:

- X Some data has been manually entered or modified. Transcription errors are possible.
- P One or more results are 'preliminary'. The data is not ready for final reporting.
- 2 There were two or more results for one analyte on one planchet imported at one time. The results in DVD may not be the same as on the raw data sheets.

Other lab defined qualifiers may occur. In general, these should be addressed in the SDG narrative.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 14

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

WORK SUMMARY

The Work Summary Report shows all samples, including QC samples, and all relevant analyses in one Sample Delivery Group (SDG). This report is often useful as supporting documentation for an invoice.

The following notes apply to this report:

- \* TEST is a code for the method used to measure associated analytes. Results and related information for each analyte are on the Data Sheet Report. In special cases, a test code used in the summary data section is not the same as in associated raw data. In this case, both codes are shown on the Work Summary.
- \* SUFFIX is the lab's code to distinguish multiple analyses (recounts, reworks, reanalyses) of a fraction of the sample. The suffix indicates which result is being reported. An empty suffix normally identifies the first attempt to analyze the sample.
- \* The LAB SAMPLE ID, TEST and SUFFIX uniquely identify all supporting data for a result. The Method Summary Report for each TEST has method performance data, such as yield, for each lab sample id and suffix and procedures used in the method.
- \* PLANCHET is an alternate lab identifier for work done for one test. It, combined with the TEST and SUFFIX, may be the best link to raw data.
- \* For QC samples, only analyses that directly QC some regular sample are shown. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.
- \* The SAS (Special Analytical Services) Number is a client or lab assigned code that reflects special processing for samples, such as rapid turn around. Counts of tests done are lists by SAS number since it is likely to affect prices.

REPORT GUIDES

Page 3

SUMMARY DATA SECTION

Page 15

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

DATA SHEET

The Data Sheet Report shows all results and primary supporting information for one client sample or Method Blank. This report corresponds to both the CLP Inorganics and Organics Data Sheet.

The following notes apply to this report:

- \* TEST is a code for the method used to measure an analyte. If the TEST is empty, no data is available; the analyte was not analyzed for.
- \* The LAB SAMPLE ID and TEST uniquely identify work within the Summary Data Section of a Data Package. The Work Summary and Method Summary Reports further identify raw data that underlies this work.

The Method Summary Report for each TEST has method performance data, such as yield, for each Lab Sample ID and a list of procedures used in the method.

- \* ERRORS can be labeled TOTAL or COUNT. TOTAL implies a preparation (non-counting method) error has been added, as square root of sum of squares, to the counting error denoted by COUNT. The preparation errors, which may vary by preparation batch, are shown on the Method Summary Report.
- \* A RESULT can be 'N.R.' (Not Reported). This means the lab did this work but chooses not to report it now, possibly because it was reported at another time.
- \* When reporting a Method Blank, a RESULT can be 'N.A.' (Not Applicable). This means there is no reported client sample work in the same preparation batch as the Blank's result. This is likely to occur when the Method Blank is associated with reanalyses of selected work for a few samples in the SDG.

The following qualifiers are defined by the DVD system:

- U The RESULT is less than the MDA (Minimum Detectable Activity).

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

REPORT GUIDES

Page 4

SUMMARY DATA SECTION

Page 16

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

DATA SHEET

If the MDA is blank, the ERROR is used as the limit.

J The RESULT is less than the RDL (Required Detection Limit) and no U qualifier is assigned.

B A Method Blank associated with this sample had a result without a U flag and, after correcting for possibly different aliquots, that result is greater than or equal to the MDA for this sample.

Normally, B is not assigned if U is. When method blank subtraction is shown on this report, B flags are assigned based on the unsubtracted values while U's are assigned based on the subtracted ones. Both flags can be assigned in this case.

For each sample result, all Method Blank results in the same preparation batch are compared. The Method Summary Report documents this and other QC relationships.

L Some Lab Control Sample that QC's this sample had a low recovery. The lab can disable assignment of this qualifier.

H Similar to 'L' except the recovery was high.

P The RESULT is 'preliminary'.

X Some data necessary to compute the RESULT, ERROR or MDA was manually entered or modified.

2 There were two or more results available for this analyte. The reported result may not be the same as in the raw data.

Other qualifiers are lab defined. Definitions should be in the SDG narrative.

The following values are underlined to indicate possible problems:

\* An MDA is underlined if it is bigger than its RDL.

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

DATA SHEET

- \* An ERROR is underlined if the 1.645 sigma counting error is bigger than both the MDA and the RESULT, implying that the MDA may not be a good estimate of the 'real' minimum detectable activity.
- \* A negative RESULT is underlined if it is less than the negative of its 2 sigma counting ERROR.
- \* When reporting a Method Blank, a RESULT is underlined if greater than its MDA. If the MDA is blank, the 2 sigma counting error is used in the comparison.

REPORT GUIDES

Page 6

SUMMARY DATA SECTION

Page 18

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

LAB CONTROL SAMPLE

The Lab Control Sample Report shows all results, recoveries and primary supporting information for one Lab Control Sample.

The following notes apply to this report:

- \* All fields in common with the Data Sheet Report have similar usage. Refer to its Report Guide for details.
- \* An amount ADDED is the lab's value for the actual amount spiked into this sample with its ERROR an estimate of the error of this amount.

An amount added is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- \* REC (Recovery) is RESULT divided by ADDED expressed as a percent.
- \* The first, computed limits for the recovery reflect:
  1. The error of RESULT, including that introduced by rounding the result prior to printing.If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.
- 2. The error of ADDED.
- 3. A lab specified, per analyte bias. The bias changes the center of the computed limits.

- \* The second limits are protocol defined upper and lower QC limits for the recovery.
- \* The recovery is underlined if it is outside either of these ranges.

REPORT GUIDES  
Page 7  
SUMMARY DATA SECTION  
Page 19

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

DUPLICATE

The Duplicate Report shows all results, differences and primary supporting information for one Duplicate and associated Original sample.

The following notes apply to this report:

- \* All fields in common with the Data Sheet Report have similar usage. This applies both to the Duplicate and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Duplicate has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- \* The RPD (Relative Percent Difference) is the absolute value of the difference of the RESULTS divided by their average expressed as a percent.

If both RESULTS are less than their MDAs, no RPD is computed and a '-' is printed.

For an analyte, if the lab did work for both samples but has data for only one, the MDA from the sample with data is used as the other's result in the RPD.

- \* The first, computed limit is the sum, as square root of sum of squares, of the errors of the results divided by the average result as a percent, hence the relative error of the difference rather than the error of the relative difference. The errors include those introduced by rounding the RESULTS prior to printing.

If this limit is labeled TOT, it includes the preparation error in the RESULTS. If labeled CNT, it does not.

This value reported for this limit is at most 999.

- \* The second limit for the RPD is the larger of:
  1. A fixed percentage specified in the protocol.

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

DUPLICATE

2. A protocol factor (typically 2) times the average MDA as a percent of the average result. This limit applies when the results are close to the MDAs.

- \* The RPD is underlined if it is greater than either limit.
- \* If specified by the lab, the second limit column is replaced by the Difference Error Ratio (DER), which is the absolute value of the difference of the results divided by the quadratic sum of their one sigma errors, the same errors as used in the first limit.

Except for differences due to rounding, the DER is the same as the RPD divided by the first RPD limit with the limit scaled to 1 sigma.

- \* The DER is underlined if it is greater than the sigma factor, typically 2 or 3, shown in the header for the first RPD limit.

REPORT GUIDES

Page 9

SUMMARY DATA SECTION

Page 21

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

MATRIX SPIKE

The Matrix Spike Report shows all results, recoveries and primary supporting information for one Matrix Spike and associated Original sample.

The following notes apply to this report:

- \* All fields in common with the Data Sheet Report have similar usage. This applies both to the Spiked and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Spike has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- \* An amount ADDED is the lab's value for the actual amount spiked into the Spike sample with its ERROR an estimate of the error of this amount.

An amount is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- \* REC (Recovery) is the Spike RESULT minus the Original RESULT divided by ADDED expressed as a percent.

- \* The first, computed limits for the recovery reflect:

1. The errors of the two RESULTS, including those introduced by rounding them prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.

2. The error of ADDED.

3. A lab specified, per analyte bias. The bias changes the center of the computed limits.

- \* The second limits are protocol defined upper and lower QC limits

REPORT GUIDES

Page 10

SUMMARY DATA SECTION

Page 22

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

MATRIX SPIKE

for the recovery.

These limits are left blank if the Original RESULT is more than a protocol defined factor (typically 4) times ADDED. This is a way of accounting for that when the spike is small compared to the amount in the original sample, the recovery is unreliable.

\* The recovery is underlined (out of spec) if it is outside either of these ranges.

REPORT GUIDES

Page 11

SUMMARY DATA SECTION

Page 23

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
Contact N. Joseph Verville

REPORT GUIDE

Client CHPRC  
Contract No. 33677  
Case no SDG H4384

METHOD SUMMARY

The Method Summary Report has two tables. One shows up to five results measured using one method. The other has performance data for the method. There is one report for each TEST, as used on the Data Sheet Report.

The following notes apply to this report:

- \* Each table is subdivided into sections, one for each preparation batch. A preparation batch is a group of aliquots prepared at roughly the same time in one work area of the lab using the same method.

There should be Lab Control Sample and Method Blank results in each preparation batch since this close correspondence makes the QC meaningful. Depending on lab policy, Duplicates need not occur in each batch since they QC sample dependencies such as matrix effects.

- \* The RAW TEST column shows the test code used in the raw data to identify a particular analysis if it is different than the test code in the header of the report. This occurs in special cases due to method specific details about how the lab labels work.

The Lab Sample or Planchet ID combined with the (Raw) Test Code and Suffix uniquely identify the raw data for each analysis.

- \* If a result is less than both its MDA and RDL, it is replaced by just 'U' on this report. If it is greater than or equal to the RDL but less than the MDA, the result is shown with a 'U' flag.

The J and X flags are as on the data sheet.

- \* Non-U results for Method Blanks are underlined to indicate possible contamination of other samples in the preparation batch. The Method Blank Report has supporting data.
- \* Lab Control Sample and Matrix Spike results are shown as: ok, No data, LOW or HIGH, with the last two underlined. 'No data'

Lab id EBRLNE  
Protocol CHPRC  
Version Ver 1.0  
Form DVD-RG  
Version 3.06  
Report date 10/04/10

REPORT GUIDES

Page 12

SUMMARY DATA SECTION

Page 24

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

METHOD SUMMARY

means no amount ADDED was specified. 'LOW' and 'HIGH' correspond to when the recovery is underlined on the Lab Control Sample or Matrix Spike Report. See these reports for supporting data.

- \* Duplicate sample results are shown as: ok, No data, or OUT, with the last two underlined. 'No data' means there was no original sample data found for this duplicate. 'OUT' corresponds to when the RPD is underlined on the Duplicate Report. See this report for supporting data.
- \* If the MDA column is labeled 'MAX MDA', there was more than one result measured by the reported method and the MDA shown is the largest MDA. If not all these results have the same RDL, the MAX MDA reflects only those results with RDL equal to the smallest one.

MDAs are underlined if greater than the printed RDL.

- \* Aliquots are underlined if less than the nominal value specified for the method.
- \* Preparation factors are underlined if greater than the nominal value specified for the method.
- \* Dilution factors are underlined if greater than the nominal value specified for the method.
- \* Residues are underlined if outside the range specified for the method. Residues are not printed if yields are.
- \* Yields, which may be gravimetric, radiometric or some type of recovery depending on the method, are underlined if outside the range specified for the method.
- \* Efficiencies are underlined if outside the range specified for the method. Efficiencies are detector and geometry dependent so this test is only approximate.

REPORT GUIDES

Page 13

SUMMARY DATA SECTION

Page 25

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.05  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899  
 Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC  
 Contract No. 33677  
 Case no SDG H4384

METHOD SUMMARY

- \* Count times are underlined if less than the nominal value specified for the method.
- \* Resolutions (as FWHM; Full Width at Half Max) are underlined if greater than the method specified limit.
- \* Tracer drifts are underlined if their absolute values are greater than the method specified limit. Tracer drifts are not printed if percent moistures are.
- \* Days Held are underlined if greater than the holding time specified in the protocol.
- \* Analysis dates are underlined if before their planchet's preparation date or, if a limit is specified, too far after it.

For some methods, ratios as percentages and error estimates for them are computed for pairs of results. A ratio column header like '1+3' means the ratio of the first result column and the third result column.

Ratios are not computed for Lab Control Sample, Method Blank or Matrix Spike results since their matrices are not necessarily similar to client samples'.

The error estimate for a ratio of results from one planchet reflects only counting errors since other errors should be correlated. For a ratio involving different planchets, if QC limits are computed based on total errors, the error for the ratio allows for the preparation errors for the planchets.

The ratio is underlined (out of spec) if the absolute value of its difference from the nominal value is greater than its error estimate. If no nominal value is specified, this test is not done.

For Gross Alpha or Gross Beta results, there may be a column showing the sum of other Alpha or Beta emitters. This sum includes all relevant

REPORT GUIDES

Page 14

SUMMARY DATA SECTION

Page 26

Lab id EBRLNE  
 Protocol CHPRC  
 Version Ver 1.0  
 Form DVD-RG  
 Version 3.06  
 Report date 10/04/10

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP H4384

SDG 7899

Contact N. Joseph Verville

GUIDE, cont.

Client CHPRC

Contract No. 33677

Case no SDG H4384

METHOD SUMMARY

results in the DVD database, whether reported or not. Results in the sum are weighted by a particles/decay value specified by the lab for each relevant analyte. Results less than their MDA are not included. No sums are computed for Lab Control, Method Blank or Matrix Spike samples since their various planchets may not be physically related.

If a ratio of total isotopic to Gross Alpha or Beta is shown, the error for the ratio reflects both the error in the Gross result and the sum, as square root of sum of squares, of the errors in the isotopic results.

For total elemental uranium or thorium results, there may be a column showing the total weight computed from associated isotopic results. Ignoring results less than their MDAs, this is a weighted sum of the isotopic results. The weights depend on the molecular weight and half-life of each isotope so as to convert activities (decays) to weight (atoms).

If a ratio of total computed to measured elemental uranium or thorium is shown, the error for the ratio reflects the errors in all the measurements.

REPORT GUIDES

Page 15

SUMMARY DATA SECTION

Page 27

Lab id EBRLNE

Protocol CHPRC

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 10/04/10

CH2M HILL Pisteau Remediation Company		CHAIN OF CUSTODY/SAMPLE ANALYSIS REQUEST		F10-275-005	PAGE 1 OF 1
COLLECTOR <i>Chamberlain</i>	COMPANY CONTACT WIDRIG, DL H4384	TELEPHONE NO. 376-2858	PROJECT COORDINATOR WIDRIG, DL	PRICE CODE 7C	DATA TURNAROUND 15 Days / 15 Days
SAMPLING LOCATION 216-N-6 EB 1 Day 1	PROJECT DESIGNATION AREA 216-N-6 Sampling and Analysis - Verification - QC Sampling	SAF NO. F10-275	COA 302222ES10	AIR QUALITY <input type="checkbox"/>	METHOD OF SHIPMENT FEDERAL EXPRESS
ICE CHEST NO. 6WS-184	FIELD LOGBOOK NO. HNF-N-507-14	ACTUAL SAMPLE DEPTH N/A	BILL OF LADING/AIR BILL NO. 7962	3788 8755	
SHIPPED TO Eberline Services	OFFSITE PROPERTY NO. 01984	SEE PTR			

MATRIX*	PRESERVATION	HMO3 to pH
A=Air		<2
DL=Drum	HOLDING TIME	6 Months
Liquids	TYPE OF CONTAINER	GP
DS=Drum	NO. OF CONTAINER(S)	1
Solids	VOLUME	1000mL
L=Liquid	SAMPLE ANALYSIS	HMO3-S3 (NCHM-63)
O=Oil	SAMPLE DATE	7-10-10
S=Soil	SAMPLE TIME	0755
SE=Sediment		✓
T=Tissue		
V=Vegetation		
W=Water		
W1=Wipe		
X=Other		

SAMPLE NO.	MATRIX*
B28157	WATER

CHAIN OF POSSESSION	SIGN/PRINT NAMES	SPECIAL INSTRUCTIONS
RELINQUISHED BY/REMOVED FROM <i>Chip Pro</i>	RECEIVED BY/STORED IN <i>V. Hall</i>	** The 200 Area S&GRP Characterization and Monitoring Sampling and Analysis GKI applies to this SAF. <input type="checkbox"/> * The laboratory shall report only the requested constituents of concern in all final data reports, and EDDs. QC results for non-requested constituents
DATE/TIME 9-10-10 1215	DATE/TIME 9-10-10 1215	
RELINQUISHED BY/REMOVED FROM <i>V. Hall</i>	RECEIVED BY/STORED IN <i>SSU</i>	
DATE/TIME 9-10-10	DATE/TIME 9-10-10	
RELINQUISHED BY/REMOVED FROM <i>SSU</i>	RECEIVED BY/STORED IN <i>LD. Wall</i>	
DATE/TIME SEP 13 2010 0800	DATE/TIME SEP 13 2010 0800	
RELINQUISHED BY/REMOVED FROM <i>Chip Pro</i>	RECEIVED BY/STORED IN <i>FEDEX</i>	
DATE/TIME SEP 13 2010 1400	DATE/TIME SEP 13 2010 1400	
RELINQUISHED BY/REMOVED FROM <i>Chip Pro</i>	RECEIVED BY/STORED IN <i>LD. Wall</i>	
DATE/TIME SEP 13 2010 1400	DATE/TIME SEP 13 2010 1400	
RELINQUISHED BY/REMOVED FROM <i>Chip Pro</i>	RECEIVED BY/STORED IN <i>LD. Wall</i>	
DATE/TIME SEP 13 2010 1400	DATE/TIME SEP 13 2010 1400	

ORIGINAL

LABORATORY SECTION	RECEIVED BY	TITLE	DATE/TIME
FINAL SAMPLE DISPOSITION	DISPOSAL METHOD	DISPOSED BY	DATE/TIME

CH2M HILL Plateau Remediation Company		CHAIN OF CUSTODY/SAMPLE ANALYSIS REQUEST		F10-275-007	PAGE 1 OF 1
COLLECTOR	FM Hill CHPRC	COMPANY CONTACT	WIDRIG, DL	PROJECT COORDINATOR	WIDRIG, DL
SAMPLING LOCATION	216-N-6 FB	TELEPHONE NO.	376-2858	SAF NO.	F10-275
ICE CHEST NO.	GWS-184	PROJECT DESIGNATION	ARRA 216-N-6 Sampling and Analysis - Verification - QC Sampling	COA	302222ES10
SHIPPED TO	Eberline Services	FIELD LOGBOOK NO.	N/A	BILL OF LADING/AIR BILL NO.	796a 3788 8755
		OFFSITE PROPERTY NO.	01984	SEE PTR	

MATRIX*	POSSIBLE SAMPLE HAZARDS/ REMARKS	PRESERVATION	HMC3 to pH
A-Air	Contains Radioactive Material at concentrations that may or may not be regulated for transportation per 49 CFR / IATA Dangerous Goods Regulations but are not releasable per DOE Order 5400.5 (1990/1993)	HOLDING TIME	<2
DL-Drum		TYPE OF CONTAINER	6 Months
DS-Drum		NO. OF CONTAINER(S)	GP
L-Liquid		VOLUME	1
O-Oil		SAMPLE ANALYSIS	1000mL
SE-Sediment			HMC4-63 (Nickel-63)
T-Tissue			
V-Vegetation			
W-Water			
WI-Wipe			
X-Other			

SAMPLE NO.	MATRIX*	SAMPLE DATE	SAMPLE TIME
B28159	WATER	9.10.10	0910

CHAIN OF POSSESSION	SIGN/ PRINT NAMES	SPECIAL INSTRUCTIONS
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	** The 200 Area S&GRP Characterization and Monitoring Sampling and Analysis GKI applies to this SAF. □* The laboratory shall report only the requested constituents of concern in all final data reports, and EDDs. QC results for non-requested constituents
FM Hill	SEP 10 2010	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
SSU-1	SEP 13 2010 0800	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 1400	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 0815	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 0815	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 0815	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 0815	
RELINQUISHED BY/REMOVED FROM	RECEIVED BY/STORED IN	
LD Wall	SEP 13 2010 0815	

ORIGINAL

LABORATORY SECTION	RECEIVED BY	TITLE	DATE/TIME
FINAL SAMPLE DISPOSITION	DISPOSAL METHOD	DISPOSED BY	DATE/TIME



**RICHMOND, CA LABORATORY**  
SAMPLE RECEIPT CHECKLIST

Client: CAPRC City METHUEN State MA  
 Date/Time received 09/14/10 0015 CoC No. FD-275--005,007  
 Container I.D. No. GWS-184 Requested TAT (Days) 15 P.O. Received Yes [ ] No [ ]

**INSPECTION**

1. Custody seals on shipping container intact? Yes  No [ ] N/A [ ]
2. Custody seals on shipping container dated & signed? Yes  No [ ] N/A [ ]
3. Custody seals on sample containers intact? Yes  No [ ] N/A [ ]
4. Custody seals on sample containers dated & signed? Yes  No [ ] N/A [ ]
5. Packing material is: Wet [ ] Dry
6. Number of samples in shipping container: 2 Sample Matrix N
7. Number of containers per sample: 1 (Or see CoC         )
8. Samples are in correct container Yes  No [ ]
9. Paperwork agrees with samples? Yes  No [ ]
10. Samples have: Tape [ ] Hazard labels [ ] Rad labels [ ] Appropriate sample labels
11. Samples are: In good condition  Leaking [ ] Broken Container [ ] Missing [ ]
12. Samples are: Preserved  Not preserved [ ] pH 1 Preservative HNO3
13. Describe any anomalies:  
\_\_\_\_\_  
\_\_\_\_\_
14. Was P.M. notified of any anomalies? Yes [ ] No [ ] Date \_\_\_\_\_
15. Inspected by [Signature] Date: 09/14/10 Time: 1230

Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	Wipe	Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	wipe
<u>fw Samples</u>	<u>260</u>						

Ion Chamber Ser. No. \_\_\_\_\_ Calibration date \_\_\_\_\_  
 Alpha Meter Ser. No. \_\_\_\_\_ Calibration date \_\_\_\_\_  
 Beta/Gamma Meter Ser. No. 99574 Calibration date 23 JUN 10