



EBERLINE SERVICES

EBERLINE ANALYTICAL CORPORATION
2030 Wright Avenue
Richmond, California 94804-3849
Phone (510) 235-2633 Fax (510) 235-0438
Toll Free (800) 841-5487
www.eberlineservices.com

May 1, 2008

Ms. Joan Kessner
Washington Closure Hanford
2620 Fermi Avenue
MSIN H4-21
Richland, WA 99352

Reference: **P.O. #S00W235A00**
Eberline Services R8-04-129-7801, SDG K1190



Dear Ms. Kessner:

Enclosed is the data report for one liquid (other liquid) sample designated under SAF No. RC-012, received at Eberline Services on April 24, 2008. The sample was analyzed according to the accompanying chain-of-custody document.

Please call if you have any questions concerning this report.

Sincerely,

For Melissa C. Mannion
Senior Program Manager

NJV

Enclosure: Data Package

1.0 GENERAL

Washington Closure Hanford (WCH) Sample Delivery Group K1190 was composed of one liquid (other liquid) sample designated under SAF No. RC-012 with a Project Designation of: 100-N Ancillary Facilities & 190-DR Waste Characterization.

The sample was received as stated on the Chain-of-Custody document. Any discrepancies are noted on the Eberline Services Sample Receipt Checklist. The results were transmitted to WCH via e-mail on May 1, 2008.

2.0 ANALYSIS NOTES

2.1 Gross Alpha and Gross Beta Analysis


No problems were encountered during the course of the analyses.

2.2 Gamma Spectroscopy

No problems were encountered during the course of the analyses.

3.0 Case Narrative Certification Statement

"I certify that this data package is in compliance with the SOW, both technically and for completeness, for other than the conditions detailed above. Release of the data obtained in this hard copy data package has been authorized by the Laboratory Manager or a designee, as verified by the following signature."

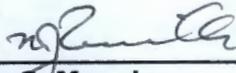

Melissa C. Mannion
Senior Program Manager

05/01/08

Date

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP K1190


SDG 7801
Contact Melissa C. Mannion

Client Hanford
Contract No. S00W235A00
Case no SDG_K1190


S U M M A R Y D A T A S E C T I O N

T A B L E O F C O N T E N T S				
About this section	.	.	.	1
Sample Summaries	.	.	.	3
Prep Batch Summary	.	.	.	5
Work Summary	.	.	.	6
Method Blanks	.	.	.	7
Lab Control Samples	.	.	.	8
Duplicates	.	.	.	9
Data Sheets	.	.	.	10
Method Summaries	.	.	.	11
Report Guides	.	.	.	14
End of Section	.	.	.	28

Prepared by

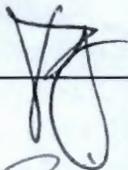

Reviewed by

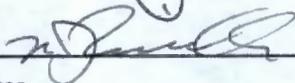

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-TOC
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. S00W235A00
Case no SDG_K1190

ABOUT THE DATA SUMMARY SECTION

The Data Summary Section of a Data Package has all data, in several useful orders, necessary for first level, routine review of the data package for a Sample Delivery Group (SDG). This section follows the Data Package Narrative, which has an overview of the data package and a discussion of special problems. It is followed by the Raw Data Section, which has full details.

The Data Summary Section has several groups of reports:

SAMPLE SUMMARIES

The Sample and QC Summary Reports show all samples, including QC samples, reported in one SDG. These reports cross-reference client and lab sample identifiers.

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches (lab groupings reflecting how work was organized) relevant to the reported SDG with information necessary to check the completeness and consistency of the SDG.

WORK SUMMARY

The Work Summary Report shows all samples and work done on them relevant to the reported SDG.

METHOD BLANKS

The Method Blank Reports, one for each Method Blank relevant to the SDG, show all results and primary supporting information for the blanks.

LAB CONTROL SAMPLES

The Lab Control Sample Reports, one for each Lab Control Sample relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 1

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. S00W235A00
Case no SDG_K1190

ABOUT THE DATA SUMMARY SECTION

DUPLICATES

The Duplicate Reports, one for each Duplicate and Original sample pair relevant to the SDG, show all results, differences and primary supporting information for these QC samples.

MATRIX SPIKES

The Matrix Spike Reports, one for each Spiked and Original sample pair relevant to the SDG, show all results, recoveries and primary supporting information for these QC samples.

DATA SHEETS

The Data Sheet Reports, one for each client sample in the SDG, show all results and primary supporting information for these samples.

METHOD SUMMARIES

The Method Summary Reports, one for each test used in the SDG, show all results, QC and method performance data for one analyte on one or two pages. (A test is a short code for the method used to do certain work to the client's specification.)

REPORT GUIDES

The Report Guides, one for each of the above groups of reports, have documentation on how to read the associated reports.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 2

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801

Contact Melissa C. Mannion

Client Hanford

Contract No. S00W235A00

Case no SDG K1190

LAB SAMPLE SUMMARY

LAB						CHAIN OF	
SAMPLE ID	CLIENT SAMPLE ID	LOCATION	MATRIX	LEVEL	SAF NO	CUSTODY	COLLECTED
R804129-01	J16HY2	1330-N Waste Pad	LIQUID		RC-012	RC-012-018	04/17/08 16:00
R804129-02	Lab Control Sample		LIQUID		RC-012		
R804129-03	Method Blank		LIQUID		RC-012		
R804129-04	Duplicate (R804129-01)	1330-N Waste Pad	LIQUID		RC-012		04/17/08 16:00

LAB SUMMARY

Page 1

SUMMARY DATA SECTION

Page 3

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-LS

Version 3.06

Report date 05/01/08

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

QC SUMMARY

Client Hanford
 Contract No. S00W235A00
 Case no SDG K1190

QC BATCH	CHAIN OF CUSTODY	CLIENT SAMPLE ID	MATRIX	% SOLIDS	SAMPLE AMOUNT	BASIS AMOUNT	DAYS SINCE RECEIVED	LAB COLL	LAB SAMPLE ID	DEPARTMENT SAMPLE ID
7801	RC-012-018	J16HY2	LIQUID		207 ml		04/24/08 7		R804129-01	7801-001
		Method Blank	LIQUID						R804129-03	7801-003
		Lab Control Sample	LIQUID						R804129-02	7801-002
		Duplicate (R804129-01)	LIQUID		207 ml		04/24/08 7		R804129-04	7801-004

QC SUMMARY

Page 1

SUMMARY DATA SECTION

Page 4

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-QS
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

PREP BATCH SUMMARY

Client Hanford
 Contract No. S00W235A00
 Case no SDG K1190

TEST	MATRIX	METHOD	PREPARATION ERROR		PLANCHETS ANALYZED				QUALI-	
			BATCH	2σ %	CLIENT	MORE	RE	BLANK		LCS
Gas Proportional Counting										
93A	LIQUID	Gross Alpha in Liquid	6148-090	20.6	1			1	1	1/1
93B	LIQUID	Gross Beta in Liquid	6148-090	11.0	1			1	1	1/1
Gamma Spectroscopy										
GAM	LIQUID	Gamma Scan in Liquid	6148-090	7.0	1			1	1	1/1

Duplicates and Matrix Spikes are those with original (Client) sample in this Sample Delivery Group.
 Blank and LCS planchets are those in the same preparation batch as some Client, Duplicate or Spike sample.

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-PBS
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

LAB WORK SUMMARY

Client Hanford
 Contract No. S00W235A00
 Case no SDG K1190

LAB SAMPLE COLLECTED RECEIVED	CLIENT SAMPLE ID LOCATION CUSTODY	MATRIX SAF No	PLANCHET	TEST	SUF- FIX	ANALYZED	REVIEWED	BY	METHOD
R804129-01 04/17/08 04/24/08	J16HY2 1330-N Waste Pad RC-012-018	LIQUID RC-012	7801-001 7801-001 7801-001	93A/93 93B/93 GAM		04/28/08 04/28/08 04/29/08	04/29/08 04/29/08 04/30/08	BW BW CSS	Gross Alpha in Liquid Gross Beta in Liquid Gamma Scan in Liquid
R804129-02	Lab Control Sample	LIQUID RC-012	7801-002 7801-002 7801-002	93A/93 93B/93 GAM		04/28/08 04/28/08 04/29/08	04/29/08 04/29/08 04/30/08	BW BW CSS	Gross Alpha in Liquid Gross Beta in Liquid Gamma Scan in Liquid
R804129-03	Method Blank	LIQUID RC-012	7801-003 7801-003 7801-003	93A/93 93B/93 GAM		04/28/08 04/28/08 04/29/08	04/29/08 04/29/08 04/30/08	BW BW CSS	Gross Alpha in Liquid Gross Beta in Liquid Gamma Scan in Liquid
R804129-04 04/17/08 04/24/08	Duplicate (R804129-01) 1330-N Waste Pad	LIQUID RC-012	7801-004 7801-004 7801-004	93A/93 93B/93 GAM		04/28/08 04/28/08 04/29/08	04/29/08 04/29/08 04/30/08	BW BW CSS	Gross Alpha in Liquid Gross Beta in Liquid Gamma Scan in Liquid

COUNTS OF TESTS BY SAMPLE TYPE

TEST	SAF No	METHOD	REFERENCE	CLIENT	MORE	RE	BLANK	LCS	DUP SPIKE	TOTAL
93A/93	RC-012	Gross Alpha in Liquid	900.0_ALPHABETA_GPC	1			1	1	1	4
93B/93	RC-012	Gross Beta in Liquid	900.0_ALPHABETA_GPC	1			1	1	1	4
GAM	RC-012	Gamma Scan in Liquid	GAMMA_GS	1			1	1	1	4
TOTALS				3			3	3	3	12

WORK SUMMARY

Page 1

SUMMARY DATA SECTION

Page 6

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-LWS
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP K1190

7801-003

Method Blank

METHOD BLANK

SDG <u>7801</u>	Client/Case no <u>Hanford</u>	SDG <u>K1190</u>
Contact <u>Melissa C. Mannion</u>	Contract <u>No. S00W235A00</u>	
Lab sample id <u>R804129-03</u>	Client sample id <u>Method Blank</u>	
Dept sample id <u>7801-003</u>	Material/Matrix <u>LIQUID</u>	
	SAF No <u>RC-012</u>	

ANALYTE	CAS NO	RESULT pCi/mL	2σ ERR (COUNT)	MDA pCi/mL	RDL pCi/mL	QUALIFIERS	TEST
Gross Alpha	12587-46-1	0.070	0.12	0.210	10.0	U	93A
Gross Beta	12587-47-2	-0.046	0.15	0.264	15.0	U	93B
Potassium 40	13966-00-2	U		0.571		U	GAM
Cobalt 60	10198-40-0	U		0.045	0.050	U	GAM
Cesium 137	10045-97-3	U		0.043	0.050	U	GAM
Radium 226	13982-63-3	U		0.099	0.100	U	GAM
Radium 228	15262-20-1	U		0.195	0.200	U	GAM
Europium 152	14683-23-9	U		<u>0.120</u>	0.100	U	GAM
Europium 154	15585-10-1	U		<u>0.143</u>	0.100	U	GAM
Europium 155	14391-16-3	U		0.088	0.100	U	GAM
Thorium 228	14274-82-9	U		0.068		U	GAM
Thorium 232	TH-232	U		0.195		U	GAM
Uranium 235	15117-96-1	U		0.163		U	GAM
Uranium 238	U-238	U		5.63		U	GAM
Americium 241	14596-10-2	U		0.053		U	GAM

100N Anclry Fclts & 190DR Wst Chrtza

QC-BLANK #65487

Lab id <u>EBRLNE</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>05/01/08</u>

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

7801-002

Lab Control Sample

LAB CONTROL SAMPLE

SDG <u>7801</u> Contact <u>Melissa C. Mannion</u> Lab sample id <u>R804129-02</u> Dept sample id <u>7801-002</u>	Client/Case no <u>Hanford</u> SDG <u>K1190</u> Contract <u>No. S00W235A00</u> Client sample id <u>Lab Control Sample</u> Material/Matrix _____ <u>LIQUID</u> SAF No <u>RC-012</u>
---	---

ANALYTE	RESULT pCi/mL	2σ ERR (COUNT)	MDA pCi/mL	RDL pCi/mL	QUALI- FIERS TEST	ADDED pCi/mL	2σ ERR pCi/mL	REC %	3σ LMTS (TOTAL)	PROTOCOL LIMITS
Gross Alpha	5.30	0.83	0.316	10.0	93A	5.10	0.20	104	59-141	70-130
Gross Beta	4.35	0.37	0.320	15.0	93B	4.67	0.19	93	80-120	70-130
Cobalt 60	2.36	0.16	<u>0.087</u>	0.050	GAM	2.19	0.088	108	83-117	80-120
Cesium 137	2.52	0.13	<u>0.095</u>	0.050	GAM	2.35	0.094	107	85-115	80-120

100N Anclyr Fclts & 190DR Wst Chrtza

QC-LCS #65486

Lab id <u>EBRLNE</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-LCS</u>
Version <u>3.06</u>
Report date <u>05/01/08</u>

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

7801-004

J16HY2

DUPLICATE

SDG <u>7801</u>	Client/Case no <u>Hanford</u>	SDG <u>K1190</u>
Contact <u>Melissa C. Mannion</u>	Contract <u>No. S00W235A00</u>	
DUPLICATE	ORIGINAL	
Lab sample id <u>R804129-04</u>	Lab sample id <u>R804129-01</u>	Client sample id <u>J16HY2</u>
Dept sample id <u>7801-004</u>	Dept sample id <u>7801-001</u>	Location/Matrix <u>1330-N Waste Pad</u> <u>LIQUID</u>
	Received <u>04/24/08</u>	Collected/Volume <u>04/17/08 16:00</u> <u>207 ml</u>
		Custody/SAF No <u>RC-012-018</u> <u>RC-012</u>

ANALYTE	DUPLICATE	2σ ERR	MDA	RDL	QUALI-	TEST	ORIGINAL	2σ ERR	MDA	QUALI-	RPD	3σ	DER
	pCi/mL	(COUNT)	pCi/mL	pCi/mL	FIERS		pCi/mL	(COUNT)	pCi/mL	FIERS	%	TOT	σ
Gross Alpha	0.034	0.083	0.148	10.0	U	93A	-0.029	0.067	0.174	U	-		1.2
Gross Beta	0.311	0.18	0.273	15.0		93B	0.087	0.16	0.262	U	113	184	1.8
Potassium 40	U		0.658		U	GAM	U		0.958	U	-		0.5
Cobalt 60	U		0.027	0.050	U	GAM	U		0.039	U	-		0.5
Cesium 137	U		0.030	0.050	U	GAM	U		0.041	U	-		0.4
Radium 226	U		0.069	0.100	U	GAM	U		0.094	U	-		0.4
Radium 228	U		0.140	0.200	U	GAM	U		<u>0.202</u>	U	-		0.5
Europium 152	U		0.096	0.100	U	GAM	U		<u>0.132</u>	U	-		0.4
Europium 154	U		0.081	0.100	U	GAM	U		<u>0.118</u>	U	-		0.5
Europium 155	U		0.089	0.100	U	GAM	U		<u>0.115</u>	U	-		0.4
Thorium 228	U		0.059		U	GAM	U		0.080	U	-		0.4
Thorium 232	U		0.140		U	GAM	U		0.202	U	-		0.5
Uranium 235	U		0.151		U	GAM	U		0.189	U	-		0.3
Uranium 238	U		3.49		U	GAM	U		5.45	U	-		0.6
Americium 241	U		0.073		U	GAM	U		0.106	U	-		0.5

100N Anclyr Felts & 190DR Wst Chrtza

QC-DUP#1 65488

DUPLICATES

Page 1

SUMMARY DATA SECTION

Page 9

Lab id <u>EBRLNE</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DUP</u>
Version <u>3.06</u>
Report date <u>05/01/08</u>

EBERLINE SERVICES / RICHMOND
SAMPLE DELIVERY GROUP K1190

7801-001

J16HY2

DATA SHEET

SDG <u>7801</u>	Client/Case no <u>Hanford</u>	<u>SDG K1190</u>
Contact <u>Melissa C. Mannion</u>	Contract No. <u>S00W235A00</u>	
Lab sample id <u>R804129-01</u>	Client sample id <u>J16HY2</u>	
Dept sample id <u>7801-001</u>	Location/Matrix <u>1330-N Waste Pad</u>	<u>LIQUID</u>
Received <u>04/24/08</u>	Collected/Volume <u>04/17/08 16:00</u>	<u>207 ml</u>
	Custody/SAF No <u>RC-012-018</u>	<u>RC-012</u>

ANALYTE	CAS NO	RESULT pCi/mL	2σ ERR (COUNT)	MDA pCi/mL	RDL pCi/mL	QUALI- FIERS	TEST
Gross Alpha	12587-46-1	-0.029	0.067	0.174	10.0	U	93A
Gross Beta	12587-47-2	0.087	0.16	0.262	15.0	U	93B
Potassium 40	13966-00-2	U		0.958		U	GAM
Cobalt 60	10198-40-0	U		0.039	0.050	U	GAM
Cesium 137	10045-97-3	U		0.041	0.050	U	GAM
Radium 226	13982-63-3	U		0.094	0.100	U	GAM
Radium 228	15262-20-1	U		<u>0.202</u>	0.200	U	GAM
Europium 152	14683-23-9	U		<u>0.132</u>	0.100	U	GAM
Europium 154	15585-10-1	U		<u>0.118</u>	0.100	U	GAM
Europium 155	14391-16-3	U		<u>0.115</u>	0.100	U	GAM
Thorium 228	14274-82-9	U		0.080		U	GAM
Thorium 232	TH-232	U		0.202		U	GAM
Uranium 235	15117-96-1	U		0.189		U	GAM
Uranium 238	U-238	U		5.45		U	GAM
Americium 241	14596-10-2	U		0.106		U	GAM

100N Anclry Fclts & 190DR Wst Chrtza

DATA SHEETS

Page 1

SUMMARY DATA SECTION

Page 10

Lab id <u>EBRLNE</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-DS</u>
Version <u>3.06</u>
Report date <u>05/01/08</u>

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

LAB METHOD SUMMARY

GROSS ALPHA IN LIQUID

GAS PROPORTIONAL COUNTING

Test 93A Matrix LIQUID

SDG 7801

Contact Melissa C. Mannion

Client Hanford

Contract No. S00W235A00

Contract SDG K1190

RESULTS

LAB	RAW	SUF-		
SAMPLE ID	TEST FIX	PLANCHET	CLIENT SAMPLE ID	Gross Alpha

Preparation batch 6148-090

R804129-01	93	7801-001	J16HY2	U
R804129-02	93	7801-002	Lab Control Sample	ok
R804129-03	93	7801-003	Method Blank	U
R804129-04	93	7801-004	Duplicate (R804129-01)	- U

Nominal values and limits from method RDLs (pCi/mL) 10.0
 100N Anclry Fclts & 190DR Wst Chrtza

METHOD PERFORMANCE

LAB	RAW	SUF-	MDA	ALIQ	PREP	DILU-	RESID	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-	
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/mL	mL	FAC	TION	mg	%	min	keV	KeV	HELD PREPARED	YZED	DETECTOR

Preparation batch 6148-090 2σ prep error 20.6 % Reference Lab Notebook #6148, pg. 90

R804129-01	93	J16HY2	0.174	2.00			20	100				11	04/28/08	04/28	GRB-109	
R804129-02	93	Lab Control Sample	0.316	2.00			60	100					04/28/08	04/28	GRB-110	
R804129-03	93	Method Blank	0.210	2.00			60	100					04/28/08	04/28	GRB-111	
R804129-04	93	Duplicate (R804129-01)	0.148	2.00			19	100					11	04/28/08	04/28	GRB-112

Nominal values and limits from method 10.0 2.00 5-250 100 180

PROCEDURES	REFERENCE	900.0_ALPHABETA_GPC
	SPP-071	Soil Dissolution, > 1.0g Aliquot, rev 5
	SPP-120	Gross Alpha and Gross Beta in Water, rev 0

AVERAGES ± 2 SD	MDA	<u>0.212</u> ± <u>0.148</u>
FOR 4 SAMPLES	RESIDUE	<u>40</u> ± <u>47</u>

METHOD SUMMARIES

Page 1

SUMMARY DATA SECTION

Page 11

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-LMS

Version 3.06

Report date 05/01/08

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

Test 93B Matrix LIQUID

SDG 7801

Contact Melissa C. Mannion

LAB METHOD SUMMARY

GROSS BETA IN LIQUID

GAS PROPORTIONAL COUNTING

Client Hanford

Contract No. S00W235A00

Contract SDG K1190

RESULTS

LAB RAW SUF-
SAMPLE ID TEST FIX PLANCHET CLIENT SAMPLE ID Gross Beta

Preparation batch 6148-090

R804129-01	93		7801-001	J16HY2	U
R804129-02	93		7801-002	Lab Control Sample	ok
R804129-03	93		7801-003	Method Blank	U
R804129-04	93		7801-004	Duplicate (R804129-01)	ok

Nominal values and limits from method RDLs (pCi/mL) 15.0
100N Anclyr Fclts & 190DR Wst Chrtza

METHOD PERFORMANCE

LAB RAW SUF- MDA ALIQ PREP DILU- RESID EFF COUNT FWHM DRIFT DAYS ANAL-
SAMPLE ID TEST FIX CLIENT SAMPLE ID pCi/mL mL PAC TION mg % min keV KeV HELD PREPARED YZED DETECTOR

Preparation batch 6148-090 2σ prep error 11.0 % Reference Lab Notebook #6148, pg. 90

R804129-01	93		J16HY2	0.262	2.00			20	100		11	04/28/08	04/28	GRB-109
R804129-02	93		Lab Control Sample	0.320	2.00			60	100			04/28/08	04/28	GRB-110
R804129-03	93		Method Blank	0.264	2.00			60	100			04/28/08	04/28	GRB-111
R804129-04	93		Duplicate (R804129-01)	0.273	2.00			19	100		11	04/28/08	04/28	GRB-112

Nominal values and limits from method 15.0 2.00 5-250 100 180

PROCEDURES	REFERENCE	900.0_ALPHABETA_GPC
	SPP-071	Soil Dissolution, > 1.0g Aliquot, rev 5
	SPP-120	Gross Alpha and Gross Beta in Water, rev 0

AVERAGES ± 2 SD	MDA	<u>0.280</u> ± <u>0.055</u>
FOR 4 SAMPLES	RESIDUE	<u>40</u> ± <u>47</u>

METHOD SUMMARIES

Page 2

SUMMARY DATA SECTION

Page 12

Lab id	<u>EBRLNE</u>
Protocol	<u>Hanford</u>
Version	<u>Ver 1.0</u>
Form	<u>DVD-LMS</u>
Version	<u>3.06</u>
Report date	<u>05/01/08</u>

EBERLINE SERVICES/RICHMOND

SAMPLE DELIVERY GROUP K1190

Test GAM Matrix LIQUID

SDG 7801

Contact Melissa C. Mannion

Client Hanford

Contract No. S00W235A00

Contract SDG K1190

LAB METHOD SUMMARY

GAMMA SCAN IN LIQUID

GAMMA SPECTROSCOPY

RESULTS

LAB	RAW	SUF-			
SAMPLE ID	TEST FIX	PLANCHET	CLIENT SAMPLE ID	Cobalt 60	Cesium 137

Preparation batch 6148-090

R804129-01	7801-001	J16HY2		U	U
R804129-02	7801-002	Lab Control Sample		ok	ok
R804129-03	7801-003	Method Blank		U	U
R804129-04	7801-004	Duplicate (R804129-01)		- U	- U

Nominal values and limits from method	RDLs (pCi/mL)	0.050	0.050
100N Anclyr Fclts & 190DR Wst Chrtza			

METHOD PERFORMANCE

LAB	RAW	SUF-	MAX MDA	ALIQ	PREP	DILU-	YIELD	EFF	COUNT	FWHM	DRIFT	DAYS	ANAL-		
SAMPLE ID	TEST FIX	CLIENT SAMPLE ID	pCi/mL	mL	FAC	TION	%	%	min	keV	KeV	HELD	PREPARED	YZED	DETECTOR

Preparation batch 6148-090			2σ prep error 7.0 %		Reference Lab Notebook #6148, pg. 90										
R804129-01	J16HY2		0.041	100					408			12	04/28/08	04/29	JR,02,00
R804129-02	Lab Control Sample		<u>0.095</u>	100					408				04/28/08	04/29	JR,05,00
R804129-03	Method Blank		0.045	100					408				04/28/08	04/29	JR,07,00
R804129-04	Duplicate (R804129-01)		0.030	100					819			12	04/28/08	04/29	JR,02,00

Nominal values and limits from method	0.050	100	100	180
---------------------------------------	-------	-----	-----	-----

PROCEDURES	REFERENCE	GAMMA_GS
	SPP-100	Ge(Li) Preparation for Commercial Samples, rev 7

AVERAGES ± 2 SD	MDA <u>0.053</u> ± <u>0.058</u>
FOR 4 SAMPLES	YIELD _____ ± _____

METHOD SUMMARIES

Page 3

SUMMARY DATA SECTION

Page 13

Lab id <u>EBRLNE</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-LMS</u>
Version <u>3.06</u>
Report date <u>05/01/08</u>

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801

Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford

Contract No. S00W235A00

Case no SDG K1190

SAMPLE SUMMARY

The Sample and QC Summary Reports show all samples, including QC samples, reported in one Sample Delivery Group (SDG).

The Sample Summary Report fully identifies client samples and gives the corresponding lab sample identification. The QC Summary Report shows at the sample level how the lab organized the samples into batches and generated QC samples. The Preparation Batch and Method Summary Reports show this at the analysis level.

The following notes apply to these reports:

- * LAB SAMPLE ID is the lab's primary identification for a sample.
- * DEPARTMENT SAMPLE ID is an alternate lab id, for example one assigned by a radiochemistry department in a lab.
- * CLIENT SAMPLE ID is the client's primary identification for a sample. It includes any sample preparation done by the client that is necessary to identify the sample.
- * QC BATCH is a lab assigned code that groups samples to be processed and QCed together. These samples should have similar matrices.

QC BATCH is not necessarily the same as SDG, which reflects samples received and reported together.

- * All Lab Control Samples, Method Blanks, Duplicates and Matrix Spikes are shown that QC any of the samples. Due to possible reanalyses, not all results for all these QC samples may be relevant to the SDG. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 14

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches in one Sample Delivery Group (SDG) with information necessary to check the completeness and consistency of the SDG.

The following notes apply to this report:

- * The preparation batches are shown in the same order as the Method Summary Reports are printed.
- * Only analyses of planchets relevant to the SDG are included.
- * Each preparation batch should have at least one Method Blank and LCS in it to validate client sample results.
- * The QUALIFIERS shown are all qualifiers other than U, J, B, L and H that occur on any analysis in the preparation batch. The Method Summary Report has these qualifiers on a per sample basis.

These qualifiers should be reviewed as follows:

- X Some data has been manually entered or modified. Transcription errors are possible.
- P One or more results are 'preliminary'. The data is not ready for final reporting.
- 2 There were two or more results for one analyte on one planchet imported at one time. The results in DVD may not be the same as on the raw data sheets.

Other lab defined qualifiers may occur. In general, these should be addressed in the SDG narrative.

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

WORK SUMMARY

The Work Summary Report shows all samples, including QC samples, and all relevant analyses in one Sample Delivery Group (SDG). This report is often useful as supporting documentation for an invoice.

The following notes apply to this report:

- * TEST is a code for the method used to measure associated analytes. Results and related information for each analyte are on the Data Sheet Report. In special cases, a test code used in the summary data section is not the same as in associated raw data. In this case, both codes are shown on the Work Summary.
- * SUFFIX is the lab's code to distinguish multiple analyses (recounts, reworks, reanalyses) of a fraction of the sample. The suffix indicates which result is being reported. An empty suffix normally identifies the first attempt to analyze the sample.
- * The LAB SAMPLE ID, TEST and SUFFIX uniquely identify all supporting data for a result. The Method Summary Report for each TEST has method performance data, such as yield, for each lab sample id and suffix and procedures used in the method.
- * PLANCHET is an alternate lab identifier for work done for one test. It, combined with the TEST and SUFFIX, may be the best link to raw data.
- * For QC samples, only analyses that directly QC some regular sample are shown. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.
- * The SAS (Special Analytical Services) Number is a client or lab assigned code that reflects special processing for samples, such as rapid turn around. Counts of tests done are lists by SAS number since it is likely to affect prices.

REPORT GUIDES

Page 3

SUMMARY DATA SECTION

Page 16

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801

Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford

Contract No. S00W235A00

Case no SDG K1190

DATA SHEET

The Data Sheet Report shows all results and primary supporting information for one client sample or Method Blank. This report corresponds to both the CLP Inorganics and Organics Data Sheet.

The following notes apply to this report:

- * TEST is a code for the method used to measure an analyte. If the TEST is empty, no data is available; the analyte was not analyzed for.
- * The LAB SAMPLE ID and TEST uniquely identify work within the Summary Data Section of a Data Package. The Work Summary and Method Summary Reports further identify raw data that underlies this work.

The Method Summary Report for each TEST has method performance data, such as yield, for each Lab Sample ID and a list of procedures used in the method.

- * ERRORS can be labeled TOTAL or COUNT. TOTAL implies a preparation (non-counting method) error has been added, as square root of sum of squares, to the counting error denoted by COUNT. The preparation errors, which may vary by preparation batch, are shown on the Method Summary Report.
- * A RESULT can be 'N.R.' (Not Reported). This means the lab did this work but chooses not to report it now, possibly because it was reported at another time.
- * When reporting a Method Blank, a RESULT can be 'N.A.' (Not Applicable). This means there is no reported client sample work in the same preparation batch as the Blank's result. This is likely to occur when the Method Blank is associated with reanalyses of selected work for a few samples in the SDG.

The following qualifiers are defined by the DVD system:

U The RESULT is less than the MDA (Minimum Detectable Activity).

REPORT GUIDES

Page 4

SUMMARY DATA SECTION

Page 17

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

GUIDE , cont .

Client Hanford
 Contract No. S00W235A00
 Case no SDG K1190

DATA SHEET

If the MDA is blank, the ERROR is used as the limit.

- J The RESULT is less than the RDL (Required Detection Limit) and no U qualifier is assigned.
 - B A Method Blank associated with this sample had a result without a U flag and, after correcting for possibly different aliquots, that result is greater than or equal to the MDA for this sample.
- Normally, B is not assigned if U is. When method blank subtraction is shown on this report, B flags are assigned based on the unsubtracted values while U's are assigned based on the subtracted ones. Both flags can be assigned in this case.
- For each sample result, all Method Blank results in the same preparation batch are compared. The Method Summary Report documents this and other QC relationships.
- L Some Lab Control Sample that QC's this sample had a low recovery. The lab can disable assignment of this qualifier.
 - H Similar to 'L' except the recovery was high.
 - P The RESULT is 'preliminary'.
 - X Some data necessary to compute the RESULT, ERROR or MDA was manually entered or modified.
 - 2 There were two or more results available for this analyte. The reported result may not be the same as in the raw data.

Other qualifiers are lab defined. Definitions should be in the SDG narrative.

The following values are underlined to indicate possible problems:

- * An MDA is underlined if it is bigger than its RDL.

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

DATA SHEET

- * An ERROR is underlined if the 1.645 sigma counting error is bigger than both the MDA and the RESULT, implying that the MDA may not be a good estimate of the 'real' minimum detectable activity.
- * A negative RESULT is underlined if it is less than the negative of its 2 sigma counting ERROR.
- * When reporting a Method Blank, a RESULT is underlined if greater than its MDA. If the MDA is blank, the 2 sigma counting error is used in the comparison.

REPORT GUIDES

Page 6

SUMMARY DATA SECTION

Page 19

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801

Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford

Contract No. S00W235A00

Case no SDG K1190

LAB CONTROL SAMPLE

The Lab Control Sample Report shows all results, recoveries and primary supporting information for one Lab Control Sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. Refer to its Report Guide for details.
- * An amount ADDED is the lab's value for the actual amount spiked into this sample with its ERROR an estimate of the error of this amount.

An amount added is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is RESULT divided by ADDED expressed as a percent.
- * The first, computed limits for the recovery reflect:
 1. The error of RESULT, including that introduced by rounding the result prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.
 2. The error of ADDED.
 3. A lab specified, per analyte bias. The bias changes the center of the computed limits.
- * The second limits are protocol defined upper and lower QC limits for the recovery.
- * The recovery is underlined if it is outside either of these ranges.

REPORT GUIDES

Page 7

SUMMARY DATA SECTION

Page 20

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
 Contract No. S00W235A00
 Case no SDG_K1190

DUPLICATE

The Duplicate Report shows all results, differences and primary supporting information for one Duplicate and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Duplicate and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Duplicate has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * The RPD (Relative Percent Difference) is the absolute value of the difference of the RESULTS divided by their average expressed as a percent.

If both RESULTS are less than their MDAs, no RPD is computed and a '-' is printed.

For an analyte, if the lab did work for both samples but has data for only one, the MDA from the sample with data is used as the other's result in the RPD.

- * The first, computed limit is the sum, as square root of sum of squares, of the errors of the results divided by the average result as a percent, hence the relative error of the difference rather than the error of the relative difference. The errors include those introduced by rounding the RESULTS prior to printing.

If this limit is labeled TOT, it includes the preparation error in the RESULTS. If labeled CNT, it does not.

This value reported for this limit is at most 999.

- * The second limit for the RPD is the larger of:
 1. A fixed percentage specified in the protocol.

REPORT GUIDES

Page 8

SUMMARY DATA SECTION

Page 21

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

DUPLICATE

2. A protocol factor (typically 2) times the average MDA as a percent of the average result. This limit applies when the results are close to the MDAs.

- * The RPD is underlined if it is greater than either limit.
- * If specified by the lab, the second limit column is replaced by the Difference Error Ratio (DER), which is the absolute value of the difference of the results divided by the quadratic sum of their one sigma errors, the same errors as used in the first limit.

Except for differences due to rounding, the DER is the same as the RPD divided by the first RPD limit with the limit scaled to 1 sigma.

- * The DER is underlined if it is greater than the sigma factor, typically 2 or 3, shown in the header for the first RPD limit.

REPORT GUIDES

Page 9

SUMMARY DATA SECTION

Page 22

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
Contract No. S00W235A00
Case no SDG_K1190

MATRIX SPIKE

The Matrix Spike Report shows all results, recoveries and primary supporting information for one Matrix Spike and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Spiked and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Spike has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * An amount ADDED is the lab's value for the actual amount spiked into the Spike sample with its ERROR an estimate of the error of this amount.

An amount is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is the Spike RESULT minus the Original RESULT divided by ADDED expressed as a percent.

- * The first, computed limits for the recovery reflect:

1. The errors of the two RESULTS, including those introduced by rounding them prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.

2. The error of ADDED.

3. A lab specified, per analyte bias. The bias changes the center of the computed limits.

- * The second limits are protocol defined upper and lower QC limits

REPORT GUIDES

Page 10

SUMMARY DATA SECTION

Page 23

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801

Contact Melissa C. Mannion

GUIDE , cont .

Client Hanford

Contract No. S00W235A00

Case no SDG_K1190

MATRIX SPIKE

for the recovery.

These limits are left blank if the Original RESULT is more than a protocol defined factor (typically 4) times ADDED. This is a way of accounting for that when the spike is small compared to the amount in the original sample, the recovery is unreliable.

- * The recovery is underlined (out of spec) if it is outside either of these ranges.

REPORT GUIDES

Page 11

SUMMARY DATA SECTION

Page 24

Lab id EBRLNE

Protocol Hanford

Version Ver 1.0

Form DVD-RG

Version 3.06

Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

REPORT GUIDE

Client Hanford
 Contract No. S00W235A00
 Case no SDG_K1190

METHOD SUMMARY

The Method Summary Report has two tables. One shows up to five results measured using one method. The other has performance data for the method. There is one report for each TEST, as used on the Data Sheet Report.

The following notes apply to this report:

- * Each table is subdivided into sections, one for each preparation batch. A preparation batch is a group of aliquots prepared at roughly the same time in one work area of the lab using the same method.

There should be Lab Control Sample and Method Blank results in each preparation batch since this close correspondence makes the QC meaningful. Depending on lab policy, Duplicates need not occur in each batch since they QC sample dependencies such as matrix effects.

- * The RAW TEST column shows the test code used in the raw data to identify a particular analysis if it is different than the test code in the header of the report. This occurs in special cases due to method specific details about how the lab labels work.

The Lab Sample or Planchet ID combined with the (Raw) Test Code and Suffix uniquely identify the raw data for each analysis.

- * If a result is less than both its MDA and RDL, it is replaced by just 'U' on this report. If it is greater than or equal to the RDL but less than the MDA, the result is shown with a 'U' flag.

The J and X flags are as on the data sheet.

- * Non-U results for Method Blanks are underlined to indicate possible contamination of other samples in the preparation batch. The Method Blank Report has supporting data.
- * Lab Control Sample and Matrix Spike results are shown as: ok, No data, LOW or HIGH, with the last two underlined. 'No data'

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
 Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
 Contract No. S00W235A00
 Case no SDG_K1190

METHOD SUMMARY

means no amount ADDED was specified. 'LOW' and 'HIGH' correspond to when the recovery is underlined on the Lab Control Sample or Matrix Spike Report. See these reports for supporting data.

- * Duplicate sample results are shown as: ok, No data, or OUT, with the last two underlined. 'No data' means there was no original sample data found for this duplicate. 'OUT' corresponds to when the RPD is underlined on the Duplicate Report. See this report for supporting data.
- * If the MDA column is labeled 'MAX MDA', there was more than one result measured by the reported method and the MDA shown is the largest MDA. If not all these results have the same RDL, the MAX MDA reflects only those results with RDL equal to the smallest one.

MDAs are underlined if greater than the printed RDL.

- * Aliquots are underlined if less than the nominal value specified for the method.
- * Preparation factors are underlined if greater than the nominal value specified for the method.
- * Dilution factors are underlined if greater than the nominal value specified for the method.
- * Residues are underlined if outside the range specified for the method. Residues are not printed if yields are.
- * Yields, which may be gravimetric, radiometric or some type of recovery depending on the method, are underlined if outside the range specified for the method.
- * Efficiencies are underlined if outside the range specified for the method. Efficiencies are detector and geometry dependent so this test is only approximate.

REPORT GUIDES

Page 13

SUMMARY DATA SECTION

Page 26

Lab id EBRLNE
 Protocol Hanford
 Version Ver 1.0
 Form DVD-RG
 Version 3.06
 Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

METHOD SUMMARY

- * Count times are underlined if less than the nominal value specified for the method.
- * Resolutions (as FWHM; Full Width at Half Max) are underlined if greater than the method specified limit.
- * Tracer drifts are underlined if their absolute values are greater than the method specified limit. Tracer drifts are not printed if percent moistures are.
- * Days Held are underlined if greater than the holding time specified in the protocol.
- * Analysis dates are underlined if before their planchet's preparation date or, if a limit is specified, too far after it.

For some methods, ratios as percentages and error estimates for them are computed for pairs of results. A ratio column header like '1+3' means the ratio of the first result column and the third result column.

Ratios are not computed for Lab Control Sample, Method Blank or Matrix Spike results since their matrices are not necessarily similar to client samples'.

The error estimate for a ratio of results from one planchet reflects only counting errors since other errors should be correlated. For a ratio involving different planchets, if QC limits are computed based on total errors, the error for the ratio allows for the preparation errors for the planchets.

The ratio is underlined (out of spec) if the absolute value of its difference from the nominal value is greater than its error estimate. If no nominal value is specified, this test is not done.

For Gross Alpha or Gross Beta results, there may be a column showing the sum of other Alpha or Beta emitters. This sum includes all relevant

REPORT GUIDES

Page 14

SUMMARY DATA SECTION

Page 27

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

EBERLINE SERVICES / RICHMOND

SAMPLE DELIVERY GROUP K1190

SDG 7801
Contact Melissa C. Mannion

GUIDE, cont.

Client Hanford
Contract No. S00W235A00
Case no SDG K1190

METHOD SUMMARY

results in the DVD database, whether reported or not. Results in the sum are weighted by a particles/decay value specified by the lab for each relevant analyte. Results less than their MDA are not included. No sums are computed for Lab Control, Method Blank or Matrix Spike samples since their various planchets may not be physically related.

If a ratio of total isotopic to Gross Alpha or Beta is shown, the error for the ratio reflects both the error in the Gross result and the sum, as square root of sum of squares, of the errors in the isotopic results.

For total elemental uranium or thorium results, there may be a column showing the total weight computed from associated isotopic results. Ignoring results less than their MDAs, this is a weighted sum of the isotopic results. The weights depend on the molecular weight and half-life of each isotope so as to convert activities (decays) to weight (atoms).

If a ratio of total computed to measured elemental uranium or thorium is shown, the error for the ratio reflects the errors in all the measurements.

REPORT GUIDES

Page 15

SUMMARY DATA SECTION

Page 28

Lab id EBRLNE
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 05/01/08

Washington Closure Hanford		CHAIN OF CUSTODY/SAMPLE ANALYSIS REQUEST						RC-012-018	Page 2 of 2 Date 4/23/08
Collector Edmundson/Isom/Perry		Company Contact Tom Edmundson		Telephone No. 376-4058		Project Coordinator KESSNER, JH		Price Code 75 9	Data Turnaround 15 Days 7-4-17-08 TRE
Project Designation 100-N Ancillary Facilities & 190-DR Waste Characterization		Sampling Location 1330-N Waste Pad 2nd Event		K1190 (7801)		SAF No. RC-012			
Ice Chest No. SML-524		Field Logbook No. EL-1516-12		COA RD4MXX2F00		Method of Shipment Fed Ex			
Shipped To EBERLINE SERVICES / LIONVILLE		Offsite Property No. A080224				Bill of Lading/Air Bill No. See OSPC			
POSSIBLE SAMPLE HAZARDS/REMARKS Potential Rad				PRESERVATION					
Special Handling and/or Storage Cool to 4°C on 4/23/08				Cool 4C	Cool 4C	None	None	None	Cool 4C
				Type of Container	aG	aG*	G/P	G/P	G/P
				No. of Container(s)	1	1	1	1	1
				Volume	1000mL	200mL	500mL 125mL	500mL 125mL	1000mL 60mL
SAMPLE ANALYSIS				PCBs - 8082	TOX - 9020	Gross Alpha	Gross Beta	See item (1) in Special Instructions	Alcohols, Glycols, & Ketones - B015 (Ethylene glycol)
				TRE 4-17-08				*	TRE 4-17-08
Sample No.	Matrix *	Sample Date	Sample Time						
J16HY1	OTHER LIQUID	4-17-08	1550			✓	✓	✓	na 4/23/08
J16HY2	OTHER LIQUID	4-17-08	1600			✓	✓	✓	
CHAIN OF POSSESSION				SIGN/PRINT NAMES				SPECIAL INSTRUCTIONS	
Relinquished By/Removed From P.R. Edmundson JR Edmundson		Date/Time 4-17-08 1820		Received By/Stored In 1060 Butte #1C		Date/Time 4-17-08 1820		(1) Gamma Spectroscopy (TCL List) (Cesium-137, Cobalt-60, Europium-152, Europium-154, Europium-155) *Gamma Spec. container is being used as shipping screen. Re-join with this set w/ screen prior to shipment to Eberline. TRE 4-17-08 Sampler unavailable to remove samples from controlled storage. Shipper removed samples from storage location taking custody of samples for shipment to lab.	
Relinquished By/Removed From 1060/1C		Date/Time 4/23/08 0930		Received By/Stored In MSTankovich		Date/Time 4/23/08			
Relinquished By/Removed From WCA MCC		Date/Time 1500 4/23/08		Received By/Stored In MSTankovich		Date/Time Feb Ex			
Relinquished By/Removed From FED EX		Date/Time		Received By/Stored In FED EX		Date/Time 04/24/08 09:30			
Relinquished By/Removed From		Date/Time		Received By/Stored In		Date/Time			
Relinquished By/Removed From		Date/Time		Received By/Stored In		Date/Time			
LABORATORY SECTION	Received By	Title				Date/Time			
FINAL SAMPLE DISPOSITION	Disposal Method	Disposed By				Date/Time			



RICHMOND, CA LABORATORY
SAMPLE RECEIPT CHECKLIST

J/K 4/24/08

Client: W.C. HANFORD City RICHLAND State WA
 Date/Time received 04/24/08 09:30 CoC No. RC-012-017
 Container I.D. No. SML-524 Requested TAT (Days) 7 P.O. Received Yes [] No []

INSPECTION

1. Custody seals on shipping container intact? Yes No [] N/A []
2. Custody seals on shipping container dated & signed? Yes No [] N/A []
3. Custody seals on sample containers intact? Yes No [] N/A []
4. Custody seals on sample containers dated & signed? Yes No [] N/A []
5. Packing material is: Wet [] Dry
6. Number of samples in shipping container: 1 Sample Matrix OTHER LIQUID
7. Number of containers per sample: 2 (Or see CoC _____)
8. Samples are in correct container Yes No []
9. Paperwork agrees with samples? Yes No []
10. Samples have: Tape [] Hazard labels [] Rad labels [] Appropriate sample labels
11. Samples are: In good condition Leaking [] Broken Container [] Missing []
12. Samples are: Preserved [] Not preserved [] pH _____ Preservative _____
13. Describe any anomalies:

14. Was P.M. notified of any anomalies? Yes [] No [] Date _____
15. Inspected by MEW Date: 04/24/08 Time: 10 30

Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	Wipe	Customer Sample No.	Beta/Gamma cpm	Ion Chamber mR/hr	wipe
<u>J164x2</u>	<u><60</u>						

Ion Chamber Ser. No. _____ Calibration date _____
 Alpha Meter Ser. No. _____ Calibration date _____
 Beta/Gamma Meter Ser. No. 100482 Calibration date 09 MAY 07