

SAF-RC-211
100-IU-2 & 100-IU-6 Miscellaneous
Restoration Sites Near 100-F –
Soil In-Process
FINAL DATA PACKAGE


COMPLETE COPY OF DATA PACKAGE TO:

Kathy Wendt

H4-21

KW 4/14/11
INITIAL/DATE

COMMENTS:

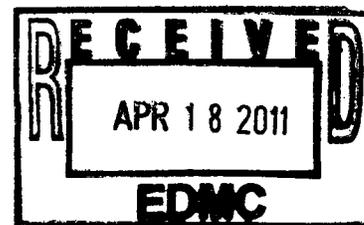

SDG D1109060

SAF-RC-211

Rad only

Chem only

Rad & Chem

Complete

Partial

Sample Location: 100-IU-2/6 600-298



ALS
Laboratory
Group
formerly
DataChem

Cover

Page 1 of 6

Report Identification Number: D1109060
 Subcontract Number: S003827A00
 Name of Industrial Hygienist: Gwen Whatley / Debbie Gothard / Ken Way
 Laboratory Identification Number: DCHM *RC-211 JK 4/1/11*
 SAF#: ~~RC-001~~ / RC-211-015
 Sample Receipt Date: 03/31/2011



Sample Information

Sample Date	Customer Sample Number	Laboratory Sample Number	Method	Analytical Batch Identification	Sample Matrix
03/29/2011	J1FXW5	1109060001	NIOSH 9002	64524	Bulk
03/29/2011	J1FXW6	1109060002	NIOSH 9002	64524	Bulk
03/29/2011	J1FXW7	1109060003	NIOSH 9002	64524	Bulk

I certify that this electronic image and all hardcopies produced from this image accurately represent the data and are in compliance with the contract specific requirements, both technically and for completeness, other than the conditions detailed above or in the sample data package narrative. Release, by submission through email, the data contained in this electronic image and the computer-readable EDD (as applicable), has been authorized by the laboratory Manager or the Manager's designee.

Name: Peter P. Steen
 Title: Chemist
 Date: April 08, 2011



Report Identification Number: D1109060
Subcontract Number: S003827A00
Name of Industrial Hygienist: Gwen Whatley / Debbie Gothard / Ken Way
Laboratory Identification Number: DCHM
SAF#: RC-001 / RC-211-015
Sample Receipt Date: 03/31/2011

General Workorder Information: There are three samples in workorder 1109059, three samples in workorder 1109060, one sample in workorder 1109165, one sample in workorder 1109167, and four samples in workorder 1109168 which were analyzed for asbestos in bulk material. No problems were encountered with the receipt of this sample.

Method Summary: All samples were examined for homogeneity. Non-homogeneous samples were ground to ensure homogeneity. Distinct layers were analyzed separately. The samples were prepared and examined for asbestos fibers utilizing the procedures outlined in NIOSH method 9002 (4th edition). A polarizing light microscope equipped with a 10x and a 16x eyepiece was used for the analysis. The area percentage of asbestos was estimated microscopically by a visual estimation of the fibers with a length-to-width aspect ratio of 3:1 or greater. If present, asbestos identities were confirmed with the appropriate refractive index oils applying dispersion staining techniques.

Sample Preparation: All samples were prepared in accordance with NIOSH method 9002 (4th edition).

Initial and Continuing Calibration Verification Analysis: N/A

Initial and Continuing Calibration Blank Analysis: N/A

Method Blank Analysis: N/A

Dilution(s): N/A.

Laboratory Control Sample and Duplicate Analysis: One Laboratory Control Sample (LCS) was prepared and analyzed with the sample batch. The results were within the control limit of +/- one reporting range.

Replicate Analysis: One sample was replicated with this analysis batch.

Flagging Codes: None

Nonconformance/Corrective Action Report (NC/CAR): N/A



ALS
Laboratory
Group
formerly
DataChem

Case Narrative

Page 3 of 6

Sample Calculation: Sample results are reported by a visual estimation of the area percentage of asbestos. If necessary, a gravimetric ashing procedure may be used to remove certain non-asbestos material from the sample; a percentage calculation is used to correct for the removal of the non-asbestos material.

Miscellaneous Comments:

- 1109059001: Brown, granular soil.
- 1109059002: Brown, granular soil.
- 1109059003: Brown, granular soil.
- 1109060001: Brown, granular soil.
- 1109060002: Brown, granular soil.
- 1109060003: Brown, granular soil.
- 1109165001: Brown, granular soil.
- 1109167001: Brown, granular soil.
- 1109168001: Brown, granular soil.
- 1109168002: Brown, granular soil.
- 1109168003: Brown, granular soil.
- 1109168004: Brown, granular soil.



ALS
Laboratory
Group
formerly
DataChem

Results

Report Identification Number: D1109060
 Subcontract Number: S003827A00
 Name of Industrial Hygienist: Gwen Whatley / Debbie Gothard / Ken Way
 Laboratory Identification Number: DCHM
 SAF#: RC-001 / RC-211-015
 Sample Receipt Date: 03/31/2011

Customer Sample Number	Laboratory Sample Number	Date Analyzed	Chrysotile %		Amosite %		Crocidolite %	
J1FXW5	1109060001	04/08/2011	<1	U	<1	U	<1	U
J1FXW6	1109060002	04/08/2011	<1	U	<1	U	<1	U
J1FXW7	1109060003	04/08/2011	<1	U	<1	U	<1	U
Required Detection Limit (RDL)			1		1		1	

Customer Sample Number	Laboratory Sample Number	Date Analyzed	Actinolite/Tremolite %		Anthophyllite %	
J1FXW5	1109060001	04/08/2011	<1	U	<1	U
J1FXW6	1109060002	04/08/2011	<1	U	<1	U
J1FXW7	1109060003	04/08/2011	<1	U	<1	U
Required Detection Limit (RDL)			1		1	

U - Parameter not detected above LOD
 J - Parameter between LOD and RDL
 ** - Not provided or unable to calculate
 NA - Not Applicable



QC Summary

Report Identification Number: D1109060
 Subcontract Number: S003827A00
 Name of Industrial Hygienist: Gwen Whatley / Debbie Gothard / Ken Way
 Laboratory Identification Number: DCHM
 SAF: RC-001 / RC-211-015
 Sample Receipt Date: 03/31/2011

Batch ID: 64524

QC Sample ID	QC Type	Analyte	Units	Result	Target
QC10788	LCS	Amosite	%	ND	ND
QC10788	LCSD	Amosite	%	ND	ND
QC10788	LCS	Chrysotile	%	20	30
QC10788	LCSD	Chrysotile	%	30	30

MB - Method Blank
 LCS - Laboratory Control Sample
 LCSD - Laboratory Control Sample Duplicate
 MS - Matrix Spike
 MSD - Matrix Spike Duplicate
 LD - Laboratory Duplicate

NA - Not Applicable
 ND - Parameter not detected above LOD

LCS, LCSD Percent Rec. = (Result / Target) * 100.0
 MS, MSD Percent Rec. = ((Result - Parent) / Target) * 100.0

LCS, LCSD Relative Percent Diff. = ((|LCS - LCSD|) / ((LCS + LCSD)/2.0)) * 100.
 MS, MSD Relative Percent Diff. = ((|MS - MSD|) / ((MS + MSD)/2.0)) * 100.
 LD Relative Percent Diff. = ((|Parent - LD|) / ((Parent + LD)/2.0)) * 100

