

0050864

Thermo Nutech
W.O. No. N9-02-056-7090

Bechtel Hanford Inc.
SDG H0346

Case Narrative

1.0 GENERAL


Bechtel Hanford Inc. Sample Delivery Group H0346 is comprised of a single liquid (water) sample designated under SAF No. C99-018 with a Project Designation of: 100-HR3-IAM (1&2) GW Monitoring Feb. 1999.

The sample was received as stated on the Chain-of-Custody document.

2.0 ANALYSIS NOTES

2.1 Total Strontium Analyses

The aliquot for the analysis was 0.5 L, less than the nominal 1.0 L however all MDA's were less than the RDL.

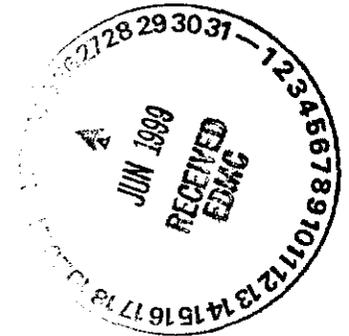

TMA/RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

SAMPLE SUMMARY

CLIENT SAMPLE ID	LOCATION	MATRIX	LEVEL	LAB SAMPLE ID	SAF NO	CHAIN OF CUSTODY	COLLECTED
B0TM35	HANFORD SITE	WATER		N902056-01	C99-018		02/12/99 08:30
Method Blank		WATER		N902056-03	C99-018		
Lab Control Sample		WATER		N902056-02	C99-018		
Duplicate (N902056-01)	HANFORD SITE	WATER		N902056-04	C99-018		02/12/99 08:30

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-CS
Version 3.06
Report date 03/22/99

TMA/RICHMOND
 SAMPLE DELIVERY GROUP H0346

SDG 7090
 Contact L.A. Johnson

Client Hanford
 Contract TRB-SBB-207925
 Case no SDG-H0346

QC SUMMARY

QC BATCH	CHAIN OF CUSTODY	CLIENT SAMPLE ID	MATRIX	SOLIDS	SAMPLE AMOUNT	BASIS AMOUNT	DAYS SINCE RECEIVED	LAB COLL SAMPLE ID	DEPARTMENT SAMPLE ID
7090		B0TM35						N902056-01	7090-001
		Method Blank						N902056-03	7090-003
		Lab Control Sample						N902056-02	7090-002
		Duplicate (N902056-01)						N902056-04	7090-004

QC SUMMARY

Page 1

SUMMARY DATA SECTION

Page 4

Lab id TMANC
 Protocol Hanford
 Version Ver 1.0
 Form DVD-QS
 Version 3.06
 Report date 03/22/99

TMA/RICHMOND

SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

PREP BATCH SUMMARY

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

TEST	MATRIX	METHOD	PREPARATION ERROR		PLANCHETS ANALYZED				QUALI-
			BATCH	2σ %	CLIENT	MORE	RE BLANK	LCS	
Beta Counting									
Y	WATER	Strontium 90 in Water	2857-159	10.0	1		1	1	1/1

Duplicates and Matrix Spikes are those with original (Client) sample in this Sample Delivery Group.
Blank and LCS planchets are those in the same preparation batch as some Client, Duplicate or Spike sample.

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-PBS
Version 3.06
Report date 03/22/99

TMA/RICHMOND

SAMPLE DELIVERY GROUP H0346

SDG 7090
 Contact L.A. Johnson

WORK SUMMARY

Client Hanford
 Contract TRB-SBB-207925
 Case no SDG-H0346

CLIENT SAMPLE ID	LAB SAMPLE ID									
LOCATION	MATRIX	COLLECTED	SUF-							
CUSTODY	SAF No	RECEIVED	PLANCHET	TEST	FIX	ANALYZED	REVIEWED	BY	METHOD	
BOTM35 HANFORD SITE	WATER	N902056-01 02/12/99	7090-001	Y		03/19/99	03/22/99	NJV	Strontium 90 in Water	
	C99-018									
Method Blank	WATER	N902056-03 02/12/99	7090-003	Y		03/19/99	03/22/99	NJV	Strontium 90 in Water	
	C99-018									
Lab Control Sample	WATER	N902056-02 02/12/99	7090-002	Y		03/19/99	03/22/99	NJV	Strontium 90 in Water	
	C99-018									
Duplicate (N902056-01) HANFORD SITE	WATER	N902056-04 02/12/99	7090-004	Y		03/19/99	03/22/99	NJV	Strontium 90 in Water	
	C99-018									

COUNTS OF TESTS BY SAMPLE TYPE

TEST	SAF No	METHOD	REFERENCE	CLIENT	MORE	RE	BLANK	LCS	DUP	SPIKE	TOTAL
Y	C99-018	Strontium 90 in Water	SR90Y90	1			1	1	1		4
TOTALS				1			1	1	1		4

Lab id TMANC
 Protocol Hanford
 Version Ver 1.0
 Form DVD-CWS
 Version 3.06
 Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

N902056-03

Method Blank

METHOD BLANK

SDG 7090 Client/Case no Hanford SDG-H0346
Contact L.A. Johnson Case no TRB-SBB-207925
Lab sample id N902056-03 Client sample id Method Blank
Dept sample id 7090-003 Material/Matrix WATER
SAF No C99-018

ANALYTE	CAS NO	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST
Strontium 90	10098-97-2	0.038	0.29	0.67	2.0	U	Y

100-HR3-IAM(1&2)GW MONITORING FEB 99

QC-BLANK 30160

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-DS
Version 3.06
Report date 03/22/99

TMA/RICHMOND

SAMPLE DELIVERY GROUP H0346

N902056-02

Lab Control Sample

LAB CONTROL SAMPLE

SDG <u>7090</u>	Client/Case no <u>Hanford</u>	<u>SDG-H0346</u>
Contact <u>L.A. Johnson</u>	Case no <u>TRB-SBB-207925</u>	
Lab sample id <u>N902056-02</u>	Client sample id <u>Lab Control Sample</u>	
Dept sample id <u>7090-002</u>	Material/Matrix <u>WATER</u>	
	SAF No <u>C99-018</u>	

ANALYTE	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS TEST	ADDED pCi/L	2σ ERR pCi/L	REC %	3σ LMTS (TOTAL)	PROTOCOL LIMITS
Strontium 90	11.9	0.39	0.39	2.0	Y	11.5	0.46	104	83-117	80-120

100-HR3-IAM(1&2)GW MONITORING FEB 99

QC-LCS 30159

Lab id <u>TMANC</u>
Protocol <u>Hanford</u>
Version <u>Ver 1.0</u>
Form <u>DVD-LCS</u>
Version <u>3.06</u>
Report date <u>03/22/99</u>

TMA/RICHMOND
SAMPLE DELIVERY GROUP H0346

N902056-04

B0TM35

DUPLICATE

SDG <u>7090</u>	Client/Case no <u>Hanford</u>	<u>SDG-H0346</u>
Contact <u>L.A. Johnson</u>	Case no <u>TRB-SBB-207925</u>	
DUPLICATE	ORIGINAL	
Lab sample id <u>N902056-04</u>	Lab sample id <u>N902056-01</u>	Client sample id <u>B0TM35</u>
Dept sample id <u>7090-004</u>	Dept sample id <u>7090-001</u>	Location/Matrix <u>HANFORD SITE</u> <u>WATER</u>
	Received _____	Collected <u>02/12/99 08:30</u>
		Custody/SAF No _____ <u>C99-018</u>

ANALYTE	DUPLICATE pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST	ORIGINAL pCi/L	2σ ERR (COUNT)	MDA pCi/L	QUALI- FIERS	RPD %	3σ TOT	PROT LIMIT
Strontium 90	26.8	1.7	1.2	2.0		Y	28.6	0.95	0.65		6	24	

100-HR3-IAM(1&2)GW MONITORING FEB 99

QC-DUP#1 30161

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-DUP
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

N902056-01

B0TM35

DATA SHEET

SDG 7090 Client/Case no Hanford SDG-H0346
Contact L.A. Johnson Case no TRB-SBB-207925
Lab sample id N902056-01 Client sample id B0TM35
Dept sample id 7090-001 Location/Matrix HANFORD SITE WATER
Received _____ Collected 02/12/99 08:30
Custody/SAF No _____ C99-018

ANALYTE	CAS NO	RESULT pCi/L	2σ ERR (COUNT)	MDA pCi/L	RDL pCi/L	QUALI- FIERS	TEST
Strontium 90	10098-97-2	28.6	0.95	0.65	2.0		Y

100-HR3-IAM(1&2)GW MONITORING FEB 99

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-DS
Version 3.06
Report date 03/22/99

TMA/RICHMOND

SAMPLE DELIVERY GROUP H0346

METHOD SUMMARY

STRONTIUM 90 IN WATER
BETA COUNTING

Test Y Matrix WATER
SDG 7090
Contact L.A. Johnson

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

RESULTS

CLIENT SAMPLE ID	LAB SAMPLE ID	RAW TEST	SUF- FIX	PLANCHET	Strontium 90
------------------	---------------	----------	----------	----------	--------------

Preparation batch 2857-159

BOTM35	N902056-01	7090-001			28.6
Method Blank	N902056-03	7090-003			U
Lab Control Sample	N902056-02	7090-002			ok
Duplicate (N902056-01)	N902056-04	7090-004			ok

Nominal values and limits from method RDLs (pCi/L) 2.0
100-HR3-IAM(1&2)GW MONITORING FEB 99

METHOD PERFORMANCE

CLIENT SAMPLE ID	LAB SAMPLE ID	RAW TEST	SUF- FIX	MDA pCi/L	ALIQ L	PREP FAC	DILU TION	YIELD %	EFF %	COUNT min	FWHM keV	DRIFT KeV	DAYS HELD	ANAL- YZED	DETECTOR
------------------	---------------	----------	----------	-----------	--------	----------	-----------	---------	-------	-----------	----------	-----------	-----------	------------	----------

Preparation batch 2857-159 2σ prep error 10.0 % Reference Lab Notebook #2857 pg. 159

BOTM35	N902056-01			0.65	<u>0.500</u>			77		200			35	03/13/99	03/19	GRB-207
Method Blank	N902056-03			0.67	1.00			73		200				03/13/99	03/19	GRB-202
Lab Control Sample	N902056-02			0.39	1.00			70		200				03/13/99	03/19	GRB-208
Duplicate (N902056-01)	N902056-04			1.2	<u>0.500</u>			81		200			35	03/13/99	03/19	GRB-203

Nominal values and limits from method 2.0 1.00 20-105 25 180

PROCEDURES	REFERENCE	SR90Y90
	EP-040	Environmental Water Dissolution, rev 1
	EP-520	Yttrium Purification for Strontium-90 Analysis, rev 0

AVERAGES ± 2 SD	MDA	<u>0.73</u> ± <u>0.68</u>
FOR 4 SAMPLES	YIELD	<u>75</u> ± <u>10</u>

METHOD SUMMARIES

Page 1

SUMMARY DATA SECTION

Page 11

Lab id	<u>TMANC</u>
Protocol	<u>Hanford</u>
Version	<u>Ver 1.0</u>
Form	<u>DVD-CMS</u>
Version	<u>3.06</u>
Report date	<u>03/22/99</u>

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

SAMPLE SUMMARY

The Sample and QC Summary Reports show all samples, including QC samples, reported in one Sample Delivery Group (SDG).

The Sample Summary Report fully identifies client samples and gives the corresponding lab sample identification. The QC Summary Report shows at the sample level how the lab organized the samples into batches and generated QC samples. The Preparation Batch and Method Summary Reports show this at the analysis level.

The following notes apply to these reports:

- * LAB SAMPLE ID is the lab's primary identification for a sample.
- * DEPARTMENT SAMPLE ID is an alternate lab id, for example one assigned by a radiochemistry department in a lab.
- * CLIENT SAMPLE ID is the client's primary identification for a sample. It includes any sample preparation done by the client that is necessary to identify the sample.
- * QC BATCH is a lab assigned code that groups samples to be processed and QCed together. These samples should have similar matrices.

QC BATCH is not necessarily the same as SDG, which reflects samples received and reported together.

- * All Lab Control Samples, Method Blanks, Duplicates and Matrix Spikes are shown that QC any of the samples. Due to possible reanalyses, not all results for all these QC samples may be relevant to the SDG. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.

REPORT GUIDES

Page 1

SUMMARY DATA SECTION

Page 12

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

PREPARATION BATCH SUMMARY

The Preparation Batch Summary Report shows all preparation batches in one Sample Delivery Group (SDG) with information necessary to check the completeness and consistency of the SDG.

The following notes apply to this report:

- * The preparation batches are shown in the same order as the Method Summary Reports are printed.
- * Only analyses of plachets relevant to the SDG are included.
- * Each preparation batch should have at least one Method Blank and LCS in it to validate client sample results.
- * The QUALIFIERS shown are all qualifiers other than U, J, B, L and H that occur on any analysis in the preparation batch. The Method Summary Report has these qualifiers on a per sample basis.

These qualifiers should be reviewed as follows:

- X Some data has been manually entered or modified. Transcription errors are possible.
- P One or more results are 'preliminary'. The data is not ready for final reporting.
- 2 There were two or more results for one analyte on one plachet imported at one time. The results in DVD may not be the same as on the raw data sheets.

Other lab defined qualifiers may occur. In general, these should be addressed in the SDG narrative.

REPORT GUIDES

Page 2

SUMMARY DATA SECTION

Page 13

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

WORK SUMMARY

The Work Summary Report shows all samples, including QC samples, and all relevant analyses in one Sample Delivery Group (SDG). This report is often useful as supporting documentation for an invoice.

The following notes apply to this report:

- * TEST is a code for the method used to measure associated analytes. Results and related information for each analyte are on the Data Sheet Report. In special cases, a test code used in the summary data section is not the same as in associated raw data. In this case, both codes are shown on the Work Summary.
- * SUFFIX is the lab's code to distinguish multiple analyses (recounts, reworks, reanalyses) of a fraction of the sample. The suffix indicates which result is being reported. An empty suffix normally identifies the first attempt to analyze the sample.
- * The LAB SAMPLE ID, TEST and SUFFIX uniquely identify all supporting data for a result. The Method Summary Report for each TEST has method performance data, such as yield, for each lab sample id and suffix and procedures used in the method.
- * PLANCHET is an alternate lab identifier for work done for one test. It, combined with the TEST and SUFFIX, may be the best link to raw data.
- * For QC samples, only analyses that directly QC some regular sample are shown. The Lab Control Sample, Method Blank, Duplicate, Matrix Spike and Method Summary Reports detail these relationships.
- * The SAS (Special Analytical Services) Number is a client or lab assigned code that reflects special processing for samples, such as rapid turn around. Counts of tests done are lists by SAS number since it is likely to affect prices.

REPORT GUIDES

Page 3

SUMMARY DATA SECTION

Page 14

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

DATA SHEET

The Data Sheet Report shows all results and primary supporting information for one client sample or Method Blank. This report corresponds to both the CLP Inorganics and Organics Data Sheet.

The following notes apply to this report:

- * TEST is a code for the method used to measure an analyte. If the TEST is empty, no data is available; the analyte was not analyzed for.
- * The LAB SAMPLE ID and TEST uniquely identify work within the Summary Data Section of a Data Package. The Work Summary and Method Summary Reports further identify raw data that underlies this work.

The Method Summary Report for each TEST has method performance data, such as yield, for each Lab Sample ID and a list of procedures used in the method.

- * ERRORS can be labeled TOTAL or COUNT. TOTAL implies a preparation (non-counting method) error has been added, as square root of sum of squares, to the counting error denoted by COUNT. The preparation errors, which may vary by preparation batch, are shown on the Method Summary Report.
- * A RESULT can be 'N.R.' (Not Reported). This means the lab did this work but chooses not to report it now, possibly because it was reported at another time.
- * When reporting a Method Blank, a RESULT can be 'N.A.' (Not Applicable). This means there is no reported client sample work in the same preparation batch as the Blank's result. This is likely to occur when the Method Blank is associated with reanalyses of selected work for a few samples in the SDG.

The following qualifiers are defined by the DVD system:

- U The RESULT is less than the MDA (Minimum Detectable Activity).

REPORT GUIDES

Page 4

SUMMARY DATA SECTION

Page 15

Lab id TMNC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

DATA SHEET

If the MDA is blank, the ERROR is used as the limit.

- J The RESULT is less than the RDL (Required Detection Limit) and no U qualifier is assigned.
- B A Method Blank associated with this sample had a result without a U flag and, after correcting for possibly different aliquots, that result is greater than or equal to the MDA for this sample.
- Normally, B is not assigned if U is. When method blank subtraction is shown on this report, B flags are assigned based on the unsubtracted values while U's are assigned based on the subtracted ones. Both flags can be assigned in this case.
- For each sample result, all Method Blank results in the same preparation batch are compared. The Method Summary Report documents this and other QC relationships.
- L Some Lab Control Sample that QC's this sample had a low recovery. The lab can disable assignment of this qualifier.
- H Similar to 'L' except the recovery was high.
- P The RESULT is 'preliminary'.
- X Some data necessary to compute the RESULT, ERROR or MDA was manually entered or modified.
- 2 There were two or more results available for this analyte. The reported result may not be the same as in the raw data.

Other qualifiers are lab defined. Definitions should be in the SDG narrative.

The following values are underlined to indicate possible problems:

- * An MDA is underlined if it is bigger than its RDL.

REPORT GUIDES

Page 5

SUMMARY DATA SECTION

Page 16

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

DATA SHEET

- * An ERROR is underlined if the 1.645 sigma counting error is bigger than both the MDA and the RESULT, implying that the MDA may not be a good estimate of the 'real' minimum detectable activity.
- * A negative RESULT is underlined if it is less than the negative of its 2 sigma counting ERROR.
- * When reporting a Method Blank, a RESULT is underlined if greater than its MDA. If the MDA is blank, the 2 sigma counting error is used in the comparison.

REPORT GUIDES

Page 6

SUMMARY DATA SECTION

Page 17

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

REPORT GUIDE

LAB CONTROL SAMPLE

The Lab Control Sample Report shows all results, recoveries and primary supporting information for one Lab Control Sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. Refer to its Report Guide for details.
- * An amount ADDED is the lab's value for the actual amount spiked into this sample with its ERROR an estimate of the error of this amount.

An amount added is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is RESULT divided by ADDED expressed as a percent.
- * The first, computed limits for the recovery reflect:
 1. The error of RESULT, including that introduced by rounding the result prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.
 2. The error of ADDED.
 3. A lab specified, per analyte bias. The bias changes the center of the computed limits.
- * The second limits are protocol defined upper and lower QC limits for the recovery.
- * The recovery is underlined if it is outside either of these ranges.

REPORT GUIDES

Page 7

SUMMARY DATA SECTION

Page 18

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

DUPLICATE

The Duplicate Report shows all results, differences and primary supporting information for one Duplicate and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Duplicate and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Duplicate has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * The RPD (Relative Percent Difference) is the absolute value of the difference of the RESULTS divided by their average expressed as a percent.

If both RESULTS are less than their MDAs, no RPD is computed and a '-' is printed.

For an analyte, if the lab did work for both samples but has data for only one, the MDA from the sample with data is used as the other's result in the RPD.

- * The first, computed limit is the sum, as square root of sum of squares, of the errors of the results divided by the average result as a percent, hence the relative error of the difference rather than the error of the relative difference. The errors include those introduced by rounding the RESULTS prior to printing.

If this limit is labeled TOT, it includes the preparation error in the RESULTS. If labeled CNT, it does not.

This value reported for this limit is at most 999.

- * The second limit for the RPD is the larger of:

1. A fixed percentage specified in the protocol.

REPORT GUIDES

Page 8

SUMMARY DATA SECTION

Page 19

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

DUPLICATE

2. A protocol factor (typically 2) times the average MDA as a percent of the average result. This limit applies when the results are close to the MDAs.

- * The RPD is underlined if it is greater than either limit.
- * If specified by the lab, the second limit column is replaced by the Difference Error Ratio (DER), which is the absolute value of the difference of the results divided by the quadratic sum of their one sigma errors, the same errors as used in the first limit.

Except for differences due to rounding, the DER is the same as the RPD divided by the first RPD limit with the limit scaled to 1 sigma.

- * The DER is underlined if it is greater than the sigma factor, typically 2 or 3, shown in the header for the first RPD limit.

REPORT GUIDES

Page 9

SUMMARY DATA SECTION

Page 20

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

REPORT GUIDE

MATRIX SPIKE

The Matrix Spike Report shows all results, recoveries and primary supporting information for one Matrix Spike and associated Original sample.

The following notes apply to this report:

- * All fields in common with the Data Sheet Report have similar usage. This applies both to the Spiked and Original sample data. Refer to the Data Sheet Report Guide for details.

If the Spike has data for a TEST and the lab did not do this test to the Original, the Original's RESULTS are underlined.

- * An amount ADDED is the lab's value for the actual amount spiked into the Spike sample with its ERROR an estimate of the error of this amount.

An amount is underlined if its ratio to the corresponding RDL is outside protocol specified limits.

- * REC (Recovery) is the Spike RESULT minus the Original RESULT divided by ADDED expressed as a percent.

- * The first, computed limits for the recovery reflect:

1. The errors of the two RESULTS, including those introduced by rounding them prior to printing.

If the limits are labeled (TOTAL), they include preparation error in the result. If labeled (COUNT), they do not.

2. The error of ADDED.

3. A lab specified, per analyte bias. The bias changes the center of the computed limits.

- * The second limits are protocol defined upper and lower QC limits

REPORT GUIDES

Page 10

SUMMARY DATA SECTION

Page 21

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

MATRIX SPIKE

for the recovery.

These limits are left blank if the Original RESULT is more than a protocol defined factor (typically 4) times ADDED. This is a way of accounting for that when the spike is small compared to the amount in the original sample, the recovery is unreliable.

- * The recovery is underlined (out of spec) if it is outside either of these ranges.

REPORT GUIDES

Page 11

SUMMARY DATA SECTION

Page 22

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

REPORT GUIDE

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

METHOD SUMMARY

The Method Summary Report has two tables. One shows up to five results measured using one method. The other has performance data for the method. There is one report for each TEST, as used on the Data Sheet Report.

The following notes apply to this report:

- * Each table is subdivided into sections, one for each preparation batch. A preparation batch is a group of aliquots prepared at roughly the same time in one work area of the lab using the same method.

There should be Lab Control Sample and Method Blank results in each preparation batch since this close correspondence makes the QC meaningful. Depending on lab policy, Duplicates need not occur in each batch since they QC sample dependencies such as matrix effects.

- * The RAW TEST column shows the test code used in the raw data to identify a particular analysis if it is different than the test code in the header of the report. This occurs in special cases due to method specific details about how the lab labels work.

The Lab Sample or Planchet ID combined with the (Raw) Test Code and Suffix uniquely identify the raw data for each analysis.

- * If a result is less than both its MDA and RDL, it is replaced by just 'U' on this report. If it is greater than or equal to the RDL but less than the MDA, the result is shown with a 'U' flag.

The J and X flags are as on the data sheet.

- * Non-U results for Method Blanks are underlined to indicate possible contamination of other samples in the preparation batch. The Method Blank Report has supporting data.
- * Lab Control Sample and Matrix Spike results are shown as: ok, No data, LOW or HIGH, with the last two underlined. 'No data'

REPORT GUIDES

Page 12

SUMMARY DATA SECTION

Page 23

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

METHOD SUMMARY

means no amount ADDED was specified. 'LOW' and 'HIGH' correspond to when the recovery is underlined on the Lab Control Sample or Matrix Spike Report. See these reports for supporting data.

- * Duplicate sample results are shown as: ok, No data, or OUT, with the last two underlined. 'No data' means there was no original sample data found for this duplicate. 'OUT' corresponds to when the RPD is underlined on the Duplicate Report. See this report for supporting data.
- * If the MDA column is labeled 'MAX MDA', there was more than one result measured by the reported method and the MDA shown is the largest MDA. If not all these results have the same RDL, the MAX MDA reflects only those results with RDL equal to the smallest one.

MDAs are underlined if greater than the printed RDL.
- * Aliquots are underlined if less than the nominal value specified for the method.
- * Preparation factors are underlined if greater than the nominal value specified for the method.
- * Dilution factors are underlined if greater than the nominal value specified for the method.
- * Residues are underlined if outside the range specified for the method. Residues are not printed if yields are.
- * Yields, which may be gravimetric, radiometric or some type of recovery depending on the method, are underlined if outside the range specified for the method.
- * Efficiencies are underlined if outside the range specified for the method. Efficiencies are detector and geometry dependent so this test is only approximate.

REPORT GUIDES

Page 13

SUMMARY DATA SECTION

Page 24

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

METHOD SUMMARY

- * Count times are underlined if less than the nominal value specified for the method.
- * Resolutions (as FWHM; Full Width at Half Max) are underlined if greater than the method specified limit.
- * Tracer drifts are underlined if their absolute values are greater than the method specified limit. Tracer drifts are not printed if percent moistures are.
- * Days Held are underlined if greater than the holding time specified in the protocol.
- * Analysis dates are underlined if before their planchet's preparation date or, if a limit is specified, too far after it.

For some methods, ratios as percentages and error estimates for them are computed for pairs of results. A ratio column header like '1+3' means the ratio of the first result column and the third result column.

Ratios are not computed for Lab Control Sample, Method Blank or Matrix Spike results since their matrices are not necessarily similar to client samples'.

The error estimate for a ratio of results from one planchet reflects only counting errors since other errors should be correlated. For a ratio involving different planchets, if QC limits are computed based on total errors, the error for the ratio allows for the preparation errors for the planchets.

The ratio is underlined (out of spec) if the absolute value of its difference from the nominal value is greater than its error estimate. If no nominal value is specified, this test is not done.

For Gross Alpha or Gross Beta results, there may be a column showing the sum of other Alpha or Beta emitters. This sum includes all relevant

REPORT GUIDES

Page 14

SUMMARY DATA SECTION

Page 25

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

TMA / RICHMOND
SAMPLE DELIVERY GROUP H0346

SDG 7090
Contact L.A. Johnson

GUIDE, cont.

Client Hanford
Contract TRB-SBB-207925
Case no SDG-H0346

METHOD SUMMARY

results in the DVD database, whether reported or not. Results in the sum are weighted by a particles/decay value specified by the lab for each relevant analyte. Results less than their MDA are not included. No sums are computed for Lab Control, Method Blank or Matrix Spike samples since their various planchets may not be physically related.

If a ratio of total isotopic to Gross Alpha or Beta is shown, the error for the ratio reflects both the error in the Gross result and the sum, as square root of sum of squares, of the errors in the isotopic results.

For total elemental uranium or thorium results, there may be a column showing the total weight computed from associated isotopic results. Ignoring results less than their MDAs, this is a weighted sum of the isotopic results. The weights depend on the molecular weight and half-life of each isotope so as to convert activities (decays) to weight (atoms).

If a ratio of total computed to measured elemental uranium or thorium is shown, the error for the ratio reflects the errors in all the measurements.

REPORT GUIDES

Page 15

SUMMARY DATA SECTION

Page 26

Lab id TMANC
Protocol Hanford
Version Ver 1.0
Form DVD-RG
Version 3.06
Report date 03/22/99

Contractor Waste Management	OFF-SITE PROPERTY CONTROL	CONTROL NO. (To be obtained from PROPERTY MANAGEMENT) W99-0-0127
--------------------------------	------------------------------	------------------------------------------------------------------------

PART I - TO BE COMPLETED BY ORIGINATOR

Department Hanford Tech Services	Section Environmental Operations	Unit Sampling Services
The following items are to be shipped from		<input checked="" type="checkbox"/> Contractor <input type="checkbox"/> Vendor
Routing <u>Fed Ex Saturday</u>		<input checked="" type="checkbox"/> Prepaid <input type="checkbox"/> Collect
Shipped to Company TMA / Western 2030 Wight Ave Address Riverside, CA 94804 City State Zip Code Country	Off-site Custodian Att: Delores Sanchez	Payroll No.
On-site Custodian		

Qty.	Property No.	Description (include Manufacture Name, Model, Serial No.)	Acquisition Cost
1	Doly	Cooler # SALL-596 C.O.C # C99-016-7 Sample # B05A135 Weight # 40 lbs Samples include base & packed on wet ice	

Classified Unclassified Shipped Under DOE Contract Shipped Under Contractor's Use Permit Contract

Necessity for the off-site use of this property

- Required for Project Work. List Project No. _____
- Business Trip
- Off-site Assignment
- Shipment to Subcontractor. List Subcontract No. _____
- Other (Please specify) Samples require analysis not available on site

CERTIFICATION OF THE RADIATION MONITORING RELEASE MUST BE SECURED THE SAME DAY THAT MATERIAL IS DELIVERED TO SHIPPING.

RM Clearance for Public Release N/A	RM Survey No. N/A	Date 2-11-99
----------------------------------------	----------------------	-----------------

Location of and Contact for Property (Name/Phone No./Bldg./Area) K.J. Yang / 372-0060 / 345 Hills		
Date Ready for Shipment 2-11-99	Cost Code to be Charged 772031101 / 108000	Approximate Date This Property will be Returned N/A
Originated By K.J. Yang	Date 2-11-99	Authorized By <i>[Signature]</i>
Property Representative Signature <i>[Signature]</i>	Date 2/12/99	Property Management Approval <i>[Signature]</i>

PART II - TO BE COMPLETED BY SHIPPING

Authorized Shipping Signature <i>[Signature]</i>	Date 2-12-99
-----------------------------------------------------	-----------------

DISTRIBUTION (AFTER FINAL SIGNATURES)

White - Property Management Yellow - Shipping Green - Accounts Payable Pink - Originator Goldenrod - Property Management

Thermo NUtech - Richmond

SAMPLE RECEIPT CHECKLIST

SAMPLE RECEIPT			
Client:	<u>PNNL</u>	Date/Time received	<u>2-16-99 9:00</u>
CoC No.	<u>699-D18-7</u>		
Container I.D. No.	<u>SML 596</u>	Requested TAT (Days)	<u>45</u> P.O. Received Yes [] No [<input checked="" type="checkbox"/>]
INSPECTION			
1.	Custody seals on shipping container intact?	Yes [<input checked="" type="checkbox"/>]	No [] N/A []
2.	Custody seals on shipping container dated & signed?	Yes [<input checked="" type="checkbox"/>]	No [] N/A []
3.	Custody seals on sample containers intact?	Yes [<input checked="" type="checkbox"/>]	No [] N/A []
4.	Custody seals on sample containers dated & signed?	Yes [<input checked="" type="checkbox"/>]	No [] N/A []
5.	Cooler Temperature: _____	Packing material is:	Wet [] Dry [<input checked="" type="checkbox"/>]
6.	Number of samples in shipping container: <u>1</u>		
7.	Number of containers per sample: <u>4</u> (Or see CoC <u>1</u>)		
8.	Paperwork agrees with samples?	Yes [<input checked="" type="checkbox"/>]	No []
9.	Samples have: Tape [] Hazard labels [] Rad labels [] Appropriate sample labels [<input checked="" type="checkbox"/>]		
10.	Samples are: In good condition [<input checked="" type="checkbox"/>] Leaking [] Broken Container [] Missing []		
11.	Describe any anomalies: <u>The samples and labels are all wet</u>		
13.	Was P.M. notified of any anomalies? Yes [] No [] Date _____		
14.	Received by <u>J.P. Corrao</u> Date: <u>2-16-99</u> Time: <u>9:00</u>		
LOGIN			
TNU W.O. No.	Group No.	Client W.O. No.	
PROGRAM MANAGER			
Sample holding times exceeded?		Yes []	No []
Client Notified: Name _____		Date/time _____	